

# **CanSAS 1-D Data Format Manual, v1.0**

Not yet licensed. Intended for unrestricted public domain use.

---

---

<b>COLLABORATORS</b>
----------------------

	<i>TITLE :</i> CanSAS 1-D Data Format Manual, v1.0		
<i>ACTION</i>	<i>NAME</i>	<i>DATE</i>	<i>SIGNATURE</i>
WRITTEN BY	Pete R. Jemian	September 26, 2009	

<b>REVISION HISTORY</b>
-------------------------

NUMBER	DATE	DESCRIPTION	NAME
v1.0	2009-09-25	converted from mediawiki to DocBook/5.0	PRJ

---

# Contents

<b>1</b>	<b>Wiki Manual</b>	<b>1</b>
1.1	Disclaimer . . . . .	1
1.2	Objective . . . . .	1
1.2.1	Status . . . . .	1
1.3	General Layout of the XML Data . . . . .	1
1.3.1	Overview . . . . .	2
1.3.1.1	Basic elements of the cansas1d/1.0 standard . . . . .	2
1.3.1.2	Required XML file header . . . . .	3
1.3.2	Rules . . . . .	3
1.3.3	Compatibility of Geometry Definitions . . . . .	5
1.3.4	Converting data into the XML format . . . . .	5
1.4	Documentation and Definitions . . . . .	5
1.4.1	XML Schema . . . . .	5
1.4.2	XML Stylesheets . . . . .	5
1.4.3	Examples and Case Studies . . . . .	6
1.4.3.1	XML layout for multiple experiments . . . . .	6
1.4.4	Foreign Elements . . . . .	7
1.4.5	Support tools for Visualization & Analysis software . . . . .	7
1.4.6	Software repositories (for canSAS1d/1.0 standard) . . . . .	8
1.5	Validation of XML against the Schema . . . . .	8
1.6	Help for XML . . . . .	8
<b>2</b>	<b>Elements of the CanSAS XML standard</b>	<b>9</b>
2.1	Required XML Header . . . . .	10
2.2	SASroot element . . . . .	13
2.3	SASentry element . . . . .	15
2.4	SASdata element . . . . .	17
2.4.1	SASdata . . . . .	18
2.4.2	Idata . . . . .	18
2.4.3	Notes . . . . .	21

---

2.5	SASsample element . . . . .	22
2.5.1	geometry . . . . .	23
2.5.2	position . . . . .	23
2.5.3	orientation . . . . .	24
2.6	SASinstrument element . . . . .	25
2.7	SASsource element . . . . .	26
2.7.1	SASsource . . . . .	26
2.7.2	beam_size . . . . .	28
2.8	SAScollimation element . . . . .	29
2.8.1	SAScollimation . . . . .	29
2.8.2	aperture . . . . .	30
2.8.3	size . . . . .	30
2.9	SASdetector element . . . . .	32
2.9.1	SASdetector . . . . .	32
2.9.2	geometry . . . . .	33
2.9.3	offset . . . . .	33
2.9.4	orientation . . . . .	34
2.9.5	beam_center . . . . .	35
2.9.6	pixel_size . . . . .	35
2.10	SASprocess element . . . . .	37
2.10.1	SASprocessnote . . . . .	38
2.11	SASnote element . . . . .	39
2.12	"{any}" element . . . . .	40
<b>3</b>	<b>Bindings and Software Support</b>	<b>41</b>
3.1	Fortran binding . . . . .	41
3.1.1	Software Development Kits . . . . .	41
3.1.2	canSAS 1-D SAS XML v1.0 support . . . . .	41
3.2	IgorPro binding . . . . .	41
3.2.1	Checkout of support code in Subversion . . . . .	42
3.2.1.1	XMLutils XOP . . . . .	42
3.2.1.2	cansasXML.ipf . . . . .	43
3.2.2	Installation . . . . .	43
3.2.3	Usage Notes . . . . .	43
3.2.4	What it does . . . . .	43
3.2.4.1	data columns . . . . .	43
3.2.4.2	metadata . . . . .	44
3.2.4.3	XML foreign namespace elements . . . . .	44
3.2.4.4	XML namespace and header . . . . .	44

---

---

3.2.4.5	XML stylesheet processing-instruction is not generated . . . . .	44
3.2.5	List of Functions . . . . .	45
3.2.6	Example test case . . . . .	46
3.2.7	Graphical User Interface . . . . .	47
3.2.7.1	Irena tool suite . . . . .	47
3.3	Java JAXB binding . . . . .	47
3.3.1	JAXB . . . . .	47
3.3.2	JAXB_cansasId_reader.java: example usage in JAVA . . . . .	48
3.3.3	example: how to retrieve I(Q) . . . . .	49
3.3.3.1	GetSASdata.java . . . . .	49
3.3.3.2	test.xml . . . . .	54
3.4	Python binding . . . . .	63
3.4.1	Comments . . . . .	63
3.4.2	gnosis.xml.objectify . . . . .	63
3.4.2.1	installation . . . . .	63
3.4.2.2	quick test in Python . . . . .	63
3.4.2.3	execution of that Python code . . . . .	64
3.4.2.4	full session output . . . . .	64
3.4.2.5	Conclusion: OK . . . . .	65
3.4.3	generateDS.py . . . . .	65
3.4.3.1	Conclusion: not ready yet . . . . .	65
3.4.4	Other suggestions . . . . .	65
<b>4</b>	<b>Other matters</b>	<b>66</b>
4.1	XML Help . . . . .	66

---

# List of Figures

1.1	block diagram of minimum elements required for cansas1d/1.0 standard . . . . .	2
1.2	definition of Q geometry for small-angle scattering . . . . .	3
1.3	definition of translation and orientation geometry as viewed from the detector towards the source . . . . .	4
2.1	The SASroot element . . . . .	10
2.2	The SASroot element . . . . .	13
2.3	The SASentry element . . . . .	15
2.4	The SASdata element . . . . .	17
2.5	Q geometry . . . . .	18
2.6	The SASsample element . . . . .	22
2.7	The SASinstrument element . . . . .	25
2.8	The SASsource element . . . . .	26
2.9	The SAScollimation element . . . . .	29
2.10	The SASdetector element . . . . .	32
2.11	The SASprocess element . . . . .	37
2.12	The SASroot element . . . . .	39

# List of Tables

1.1	Basic elements of the CanSAS 1-D standard . . . . .	2
3.1	metadata for the "cs_collagen_full.xml" case study . . . . .	44



# Preface



CanSAS 1-D Data Format Manual, v1.0

---

## Note

Provide a note here describing CanSAS and the initiative to establish a standard file format.

---

"CanSAS" stands for "Collective Action for Nomadic Small-Angle Scatterers."

This work is the initiative of the CanSAS 1D Data Formats Working Group, established at the canSAS-V workshop, NIST, Gaithersburg, Maryland, USA from October 29th to 31st 2007. It derives many of its foundations from previous works such as the SASXML format, a joint collaboration between ISIS and ILL.

# Chapter 1

## Wiki Manual

### 1.1 Disclaimer

This description is meant to inform the community how to layout the information within the XML files. However, should the information in this document and the [cansas1d/1.0 SAS XML Schema](http://svn.smallangles.net/trac/canSAS/browser/ldwg/trunk/cansas1d/1.0/SAS_XML_Schema) ([http://svn.smallangles.net/trac/canSAS/browser/ldwg/trunk/cansas1d/1.0/SAS\\_XML\\_Schema](http://svn.smallangles.net/trac/canSAS/browser/ldwg/trunk/cansas1d/1.0/SAS_XML_Schema)) differ, the XML Schema will be deemed to have the most correct description of the standard.

### 1.2 Objective

One of the first aims of the canSAS (Collective Action for Nomadic Small-Angle Scatterers) forum of users, software developers, and facility staff was to discuss better sharing of SAS data analysis software. [CanSAS](http://www.smallangles.net/canSAS) (<http://www.smallangles.net/canSAS>) identified that a significant need within the SAS community can be satisfied by a robust, self-describing, text-based, standard format to communicate reduced one-dimensional small-angle scattering data, "I(Q)", between users of our facilities. Our goal has been to define such a format that leaves the data file instantly human-readable, editable in the simplest of editors, and importable by simple text import filters in programs that need not recognise advanced structure in the file nor require advanced programming interfaces. The file should contain both the primary data of "I(Q)" and also any other descriptive information (metadata) about the sample, measurement, instrument, processing, or analysis steps.

The cansas1d/1.0 standard meets the objectives for a 1D standard, incorporating metadata about the measurement, parameters and results of processing or analysis steps. Even multiple measurements (related or unrelated) may be included within a single XML file.

#### 1.2.1 Status

Version 1.0 was tagged from the subversion repository on 2009-05-12 as no changes were committed since January 2009. Use this command to checkout the tagged release.

```
svn checkout http://svn.smallangles.net/svn/canSAS/ldwg/tags/v1.0 cansas1d-1.0
```

### 1.3 General Layout of the XML Data

The canSAS 1-D standard for reduced 1-D SAS data is implemented using XML files. A single file can contain SAS data from a single experiment or multiple experiments. All types of relevant data ("I(Q)", metadata) are described for each experiment. More details are provided below.

---

### 1.3.1 Overview

The basic elements of the cansas1d/1.0 standard are shown in the following table. After an XML header, the root element of the file is `[[cansas1d_SASroot | SASroot]]` which contains one or more `[[cansas1d_SASentry | SASentry]]` elements, each of which describes a single experiment (data set, time-slice, step in a series, new sample, etc.). Details of the `[[cansas1d_SASentry | SASentry]]` element are also shown in the next figure. Refer to the `[[cansas1d_block_diagrams | block diagrams]]` for alternative depictions. See <http://svn.smallangles.net/trac/canSAS/browser/1dwg/trunk/cansas1d.xml> for an example XML file. Examples, Case Studies, and other background information are below. More discussion can be found on the `[[1D_Data_Formats_Working_Group|1D Data Formats Working Group]]` page and its `[[Talk:1D_Data_Formats_Working_Group|discussion]]` page. Details about each specific field (XPath string, XML elements and attributes) are described on the `[[cansas1d_definition_of_terms]]` page.

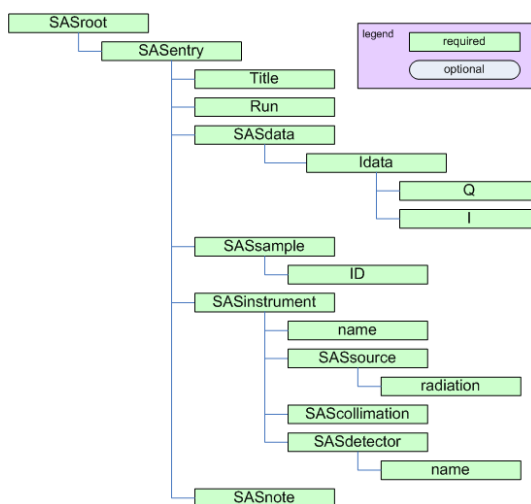


Figure 1.1: block diagram of minimum elements required for cansas1d/1.0 standard

- SASroot: the root element of the file (after the XML header)
- SASentry: describes a single experiment (data set, time-slice, step in a series, new sample, etc.)
- block diagrams
- cansas1d.xml example XML file
- discussion of this format: basic more
- Seek outside help for XML
- Definition of terms: Details about each specific field (XPath string, XML elements and attributes)

#### 1.3.1.1 Basic elements of the cansas1d/1.0 standard

Table 1.1: Basic elements of the CanSAS 1-D standard

Element	Description
XML Header	descriptive info required at the start of every XML file
SASroot	root element of XML file
SASentry	data set, time-slice, step in a series, new sample, etc.
:::Title	for this particular SASentry
:::Run	run number or ID number of experiment
:::[[cansas1d_any   {any}]]	any non-cansas1d/1.0 element can be used at this point
:::[[cansas1d_SASdata   SASdata]]	this is where the reduced 1-D SAS data is stored
:::[[cansas1d_SASdata   Idata]]	a single data point in the dataset

Table 1.1: (continued)

Element	Description
...[[cansas1d_any   {any}]]	any non-cansas1d/1.0 element can be used at this point
...[[cansas1d_SASsample   SASsample]]	description of the sample
...[[cansas1d_SASinstrument   SASinstrument]]	description of the instrument
...[[cansas1d_SASsource   SASsource]]	description of the source
...[[cansas1d_SAScollimation   SAScollimation]]	description of the collimation
...[[cansas1d_SASdetector   SASdetector]]	description of the detector
...[[cansas1d_SASprocess   SASprocess]]	for each processing or analysis step
...[[cansas1d_SASnote   SASnote]]	anything at all

### 1.3.1.2 Required XML file header

```
<?xml version="1.0"?>
<?xml-stylesheet type="text/xsl" href="cansasxml-html.xsl" ?>
<SASroot version="1.0"
xmlns="cansas1d/1.0"
xmlns:xsi="http://www.w3.org/2001/XMLSchema-instance"
xsi:schemaLocation="cansas1d/1.0
http://svn.smallangles.net/svn/canSAS/1dwg/trunk/cansas1d.xsd"
>
```

### 1.3.2 Rules

1. canSAS1d/1.0 XML data files will adhere to the standard if they can successfully [[cansas1d\_documentation#Validation\_of\_XML\_ | validate]] against the established XML Schema ([<http://svn.smallangles.net/trac/canSAS/browser/1dwg/trunk/cansas1d.xsd> cansas1d.xsd])
- 2.

$$Q = (4\pi/\lambda)\sin(\theta) \quad (1.1)$$

where  $\lambda$  is the wavelength of the radiation

and  $2\theta$  is the angle through which the detected radiation has been scattered.

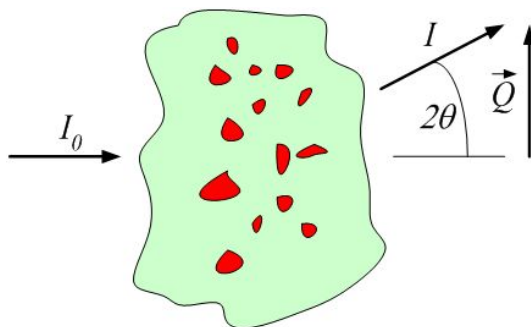


Figure 1.2: definition of Q geometry for small-angle scattering

3. units to be given in standard SI abbreviations (eg, m, cm, mm, nm, K) with the following exceptions:
  - a.  $\mu\text{m}$ =micrometres

- b. C=celsius
  - c. A=Angstroms
  - d. percent=%.
  - e. fraction
  - f. a.u.=arbitrary units
  - g. none=no units are relevant (such as dimensionless)
4. where reciprocal units need to be quoted the format shall be "1/abbreviation"
  5. when raised to a power, use similar to "A<sup>3</sup>" or "1/m<sup>4</sup>" (and not "A3" or "m-4")
  6. axes:
    - a. z is along the flight path (positive value in the direction of the detector)
    - b. x is orthogonal to z in the horizontal plane (positive values increase to the right when viewed towards the incoming radiation)
    - c. y is orthogonal to z and x in the vertical plane (positive values increase upwards)

Table 1.2:

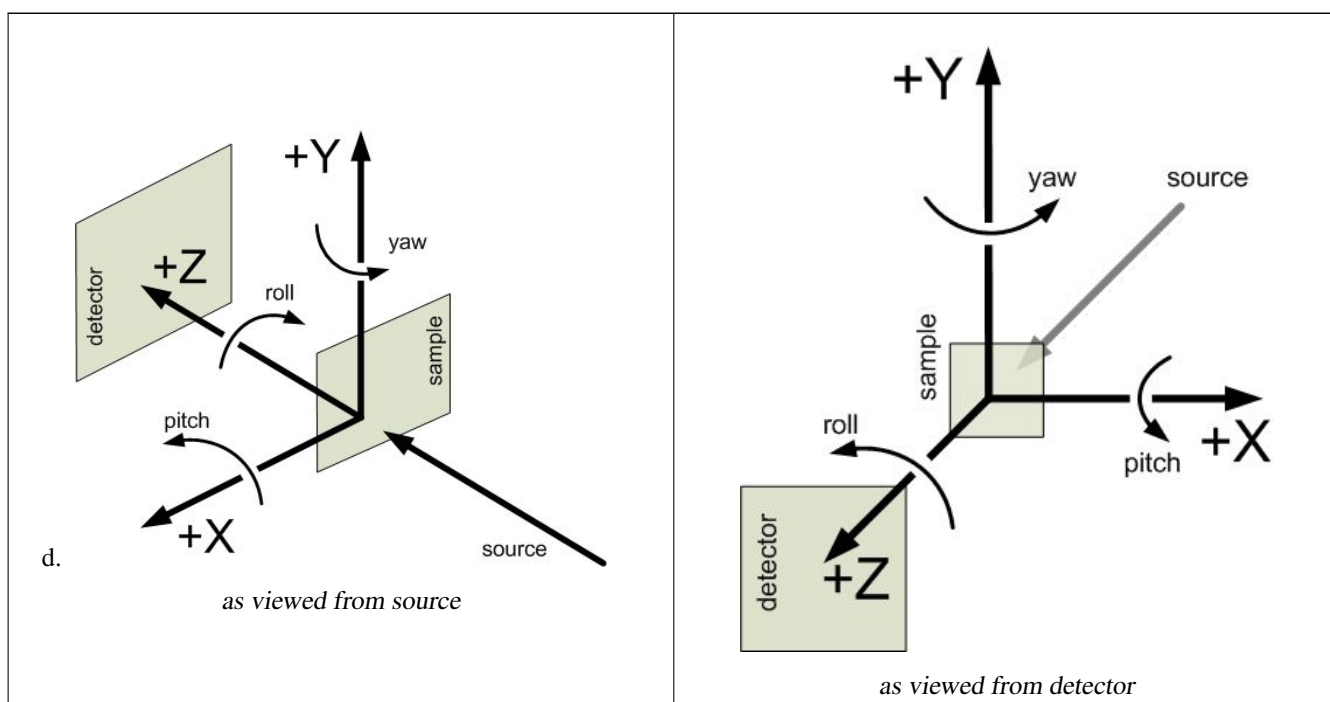


Figure 1.3: definition of translation and orientation geometry as viewed from the detector towards the source

7. orientation (angles) describes one-axis rotations (rotations about multiple axes require more information):
  - a. roll is about z
  - b. pitch is about x
  - c. yaw is about y
8. Unicode characters MUST NOT be used
9. Binary data is not supported

### 1.3.3 Compatibility of Geometry Definitions

Note: translation and orientation geometry used by canSAS are consistent with:

- [http://en.wikipedia.org/wiki/Cartesian\\_coordinate\\_system](http://en.wikipedia.org/wiki/Cartesian_coordinate_system)
- [http://en.wikipedia.org/wiki/Right-hand\\_rule](http://en.wikipedia.org/wiki/Right-hand_rule)
- [http://www.nexusformat.org/Coordinate\\_Systems](http://www.nexusformat.org/Coordinate_Systems)
- <http://mcstas.risoe.dk/documentation/tutorial/node6.html>
- <http://webhost5.nts.jhu.edu/reza/book/kinematics/kinematics.htm>

The translation and orientation geometry definitions used here are different than those used by "SHADOW" (<http://www.nanotech.wisc.edu/shadow/>) where the "y" and "z" axes are swapped and the direction of "x" is changed.

### 1.3.4 Converting data into the XML format

A WWW form (<http://www.smallangles.net/canSAS/xmlWriter/> canSAS/xmlWriter) is available to translate three-column ASCII text data into the canSAS1d/1.0 XML format. This form will help you in creating an XML file with all the required elements in the correct places. The form requests the SAS data of Q, I, and Idev (defined elsewhere on this page) and some basic metadata (title, run, sample info, ...). Press the `<nowiki>Submit</nowiki>` button and you will receive a nicely formatted WWW page with the SAS data. If you then choose "View page source" (from one of your browser menus), you will see the raw XML of the canSAS1d/1.0 XML format and you can copy/paste this into an XML file.

The SAS data that you paste into the form box is likely to be copied directly from a 3-column ASCII file from a text editor. Line breaks are OK, they will be treated as white-space as will tabs and commas. Do not be concerned that the data looks awful in the form entry box, just check the result to see that it comes out OK.

## 1.4 Documentation and Definitions

- **Documentation:** [\[\[cansas1d\\_documentation\]\]](#) (this page)
- **Definitions:** [\[\[cansas1d\\_definition\\_of\\_terms\]\]](#)
- **Block diagrams:** [\[\[cansas1d\\_block\\_diagrams\]\]](#)

### 1.4.1 XML Schema

The <http://svn.smallangles.net/trac/canSAS/browser/1dwg/trunk/cansas1d.xsd> cansas1d.xsd <http://www.w3schools.com/xsd/XMLSchema/> XML Schema] defines the rules for the XML file format (<http://svn.smallangles.net/trac/canSAS/browser/1dwg/trunk/cansas1d.xsd> TRAC], <http://svn.smallangles.net/svn/canSAS/1dwg/trunk/cansas1d.xsd> SVN]) and is used to validate any XML file for adherence to the format.

### 1.4.2 XML Stylesheets

- "cansasxml-html.xml": [http://www.w3schools.com/xsl/XSLT\\_stylesheets/](http://www.w3schools.com/xsl/XSLT_stylesheets/) can be used to extract metadata or to convert into another file format. The default canSAS stylesheet [<http://svn.smallangles.net/svn/canSAS/1dwg/trunk/cansasxml-html.xml> cansasxml-html.xml] should be copied into each folder with canSAS XML data file(s). It can be used to display the data in a supporting WWW browser (such as Firefox or Internet Explorer) or to import into Microsoft Excel (with the added XML support in Excel). (See the excellent write-up by Steve King, ISIS, at [http://www.isis.rl.ac.uk/archive/LargeScale/LOQ/xml/cansas\\_xml\\_for\\_an\\_example/](http://www.isis.rl.ac.uk/archive/LargeScale/LOQ/xml/cansas_xml_for_an_example/).) By default, MS Windows binds "\*.xml" files to start Internet Explorer. Double-clicking on a canSAS XML data file with the [<http://svn.smallangles.net/trac/canSAS/browser/1dwg/trunk/cansasxml-html.xml> "cansasxml-html.xml"] stylesheet in the same directory will produce a WWW page with the SAS data and selected metadata.

- Suggestions for support software that writes canSAS1d/1.0 XML data files:
  - be sure to update to the latest SVN repository revision (command: `svn update`)
  - check the output directory to see if it contains the default XSLT file.
  - copy the latest XSLT file to the output directory if either:
    - \* the output directory contains an older revision
    - \* the output directory does not have the default XSLT file
  - The most recent XSLT file can be identified by examining the file for the `$$$ Revision: $$$` string. For example:

```
# \$$Revision: 106 $
```

is version 66 (updated 2009-01-12).

### 1.4.3 Examples and Case Studies

- [<http://svn.smallangles.net/trac/canSAS/browser/1dwg/trunk/cansas1d.xml> cansas1d.xml] basic example: Note that, for clarity, only one row of data is shown. This is probably a very good example to use as a starting point for creating XML files with a text editor.
- [<http://svn.smallangles.net/trac/canSAS/browser/1dwg/trunk/bimodal-test1.xml> bimodal-test1.xml]: Simulated SAS data to test size distribution calculation routines.
- [[Glassy\_Carbon\_Round\_Robin | Glassy Carbon Round Robin]]: Glassy carbon samples measured at several facilities worldwide.
- [[cansas1d\_casestudy\_collagen | dry chick collagen]]: illustrates the minimum information necessary to meet the requirements of the standard format
- [[cansas1d\_casestudy\_af1410 | AF1410 steel]]: SANS study using magnetic contrast variation (with multiple samples and multiple data sets for each sample), the files can be viewed from TRAC (no description yet): [<http://svn.smallangles.net/trac/canSAS/browser/1dwg/trunk/examples/af1410/>]
- [<http://svn.smallangles.net/trac/canSAS/browser/1dwg/trunk/cansas1d-template.xml> cansas1d-template.xml]: This is used to test all the rules in the XML Schema. This is probably not a very good example to use as a starting point for creating XML files with a text editor since it tests many of the special-case rules.

#### 1.4.3.1 XML layout for multiple experiments

Each experiment is described with a single `$$$SASentry$$$` element. The fragment below shows how multiple experiments can be included in a single XML file. Full examples of canSAS XML files with multiple experiments include:

- [<http://svn.smallangles.net/trac/canSAS/browser/1dwg/trunk/W1W2.XML> ISIS LOQ SANS instrument]
- [[http://svn.smallangles.net/trac/canSAS/browser/1dwg/trunk/cs\\_af1410.xml](http://svn.smallangles.net/trac/canSAS/browser/1dwg/trunk/cs_af1410.xml) NIST SANS data]

Here is a brief sketch of how a file would be arranged with multiple `SASentry` elements and multiple `SASdata` elements.

```
<?xml version="1.0"?>
<?xml-stylesheet type="text/xsl" href="cansasxml-html.xsl" ?>
<SASroot version="1.0"
xmlns="cansas1d/1.0"
xmlns:xsi="http://www.w3.org/2001/XMLSchema-instance"
xsi:schemaLocation="cansas1d/1.0 http://svn.smallangles.net/svn/ ↵
    cansas1dwg/trunk/cansas1d.xsd"
>
<SASentry name="071121.dat#S22">
<!-- contents of the first experiment in the file go here -->
</SASentry>
```

```
<SASentry name="example temperature series">
<!-- example with two SAS data sets related to the same sample -->
</-->
<Title>title of this series</Title>
<Run name="run1">42-001</Run>
<Run name="run2">42-002</Run>
<SASdata name="run1">
<!-- data from 42-001 run comes here -->
</SASdata>
<SASdata name="run2">
<!-- data from 42-002 run comes here -->
</SASdata>
<!-- other elements come here for this entry -->
</SASentry>
<SASentry name="other sample">
<!-- any number of additional experiments can be included, as desired -->
<!-- SASentry elements in the same XML file do not have to be related -->
</SASentry>
</SASroot>
```

#### 1.4.4 Foreign Elements

To allow for inclusion of elements that are not defined by the `cansas1d.xsd` XML Schema, XML "foreign elements" are permitted at select locations in the `cansas1d/1.0` format. Please refer to the references (and others) [[#Help\_for\_XML | below]] for deeper discussions on foreign elements.

No examples exist.

---

##### Note

2009-09, PRJ: Actually, the ISIS glassy carbon data is now an example

---

At present, all examples of `canSAS` xml files using foreign namespaces have been converted to bring that data into either the `""SASprocessnote""` or `""SASnote""` elements. Refer to the [http://svn.smallangles.net/trac/canSAS/changeset/47 TRAC changes] for an example of arranging the content in `""SASprocessnote""` to avoid the use of foreign namespace elements.

#### 1.4.5 Support tools for Visualization & Analysis software

Support for importing `canSAS1d/1.0` files exists for these languages:

---

##### Note

Refactor the wiki pages here and link as appropriate.

---

- **FORTRAN**

- **IgorPro**

- **Java**

- **Microsoft Excel** Support for Microsoft Excel is provided through the default `canSAS` stylesheet [cansasxml-html.xml](#). An [excellent description](#) of how to import data from the `cansas1d/1.0` format into Excel is available from the [ISIS LOQ instrument](#) (<http://www.isis.stfc.ac.uk/instruments/loq/loq2470.html>). Also note the [old WWW site](#). (<http://www.isis.rl.ac.uk/LargeScale/LOQ/loq2470.html>) may still be available.

---



- **PHP**
- **Python**
- **XSLT** (useful in a web browser)

#### 1.4.6 Software repositories (for canSAS1d/1.0 standard)

- **TRAC:** <http://svn.smallangles.net/trac/canSAS/browser/1dwg/tags/v1.0>
- **Subversion:** <http://svn.smallangles.net/svn/canSAS/1dwg/tags/v1.0>

### 1.5 Validation of XML against the Schema

1. open browser to: <http://www.xmlvalidation.com/>
2. paste content of candidate XML file (with reference in the header to the XML Schema as shown above) into the form
3. press <validate>
4. paste content of [<http://svn.smallangles.net/svn/canSAS/1dwg/trunk/cansas1d.xsd> cansas1d.xsd] XSD file into form and press <continue validation>
5. check the results

### 1.6 Help for XML

The various references for help on XML have been moved to their own wiki page: [\[\[XmlHelp\]\]](#)

---

#### Note

chapter="Other" section="XML Help" xml:id="wiki-XML\_Help"

---

## Chapter 2

# Elements of the CanSAS XML standard

There are various elements (tag names) in the canSAS1d/1.0 standard. Each of these is described below.

2.1 Required XML Header

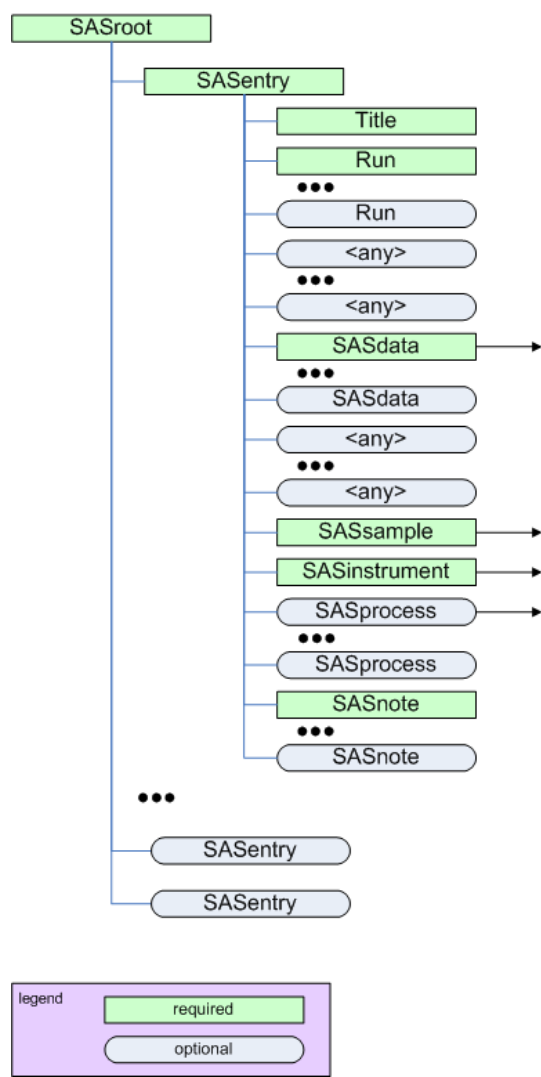


Figure 2.1: The SASroot element

Table 2.1:

Name	Type	occurrence	Description	Attributes
<i>xml declaration</i>	header	[1..1]	<?xml version="1.0"?>	version="1.0"


Table 2.1: (continued)

<i>stylesheet</i>	header	[0..1]	<pre>&lt;?xml-stylesheet type="text/xs- l" href="example- .xml" ?&gt;</pre> <p>Declares that <i>example.xml</i> (needs to be in the local directory) will be the default stylesheet to an XML visualization tool. Change <i>example.xml</i> to indicate a different stylesheet in the local directory. Refer to <a href="#">W3 Schools XSLT Help</a> for assistance in constructing XSLT files. XML rules actually allow for multiple stylesheet declarations. Explore this possibility as your own adventure.</p>	<pre>type="text/xs- l" href="example- .xml"</pre>
-------------------	--------	--------	--	---

Table 2.1: (continued)

SASroot	container	[1..1]	<p>The canSAS reduced 1-D SAS data (cansas1d/1.0) will be in the <b>SASroot</b> database. (This is similar to NXroot used by NeXus.)</p> <ul style="list-style-type: none"><li>• <b>version="1.0"</b> is required to identify the cansas1d/1.0 standard for SAS data.</li><li>• <b>xmlns</b> sets the default namespace for all elements (with no prefix) in this file.</li><li>• <b>xmlns:xsi</b> sets <b>xsi</b> as the prefix for any elements from the governing XML Schema.</li><li>• <b>xsi:schemaLocation</b> associates a suggested URL (where the cansas1d/1.0 XML Schema might be found) with the default namespace string.</li></ul> <pre>&lt;SASroot   xmlns="   xmlns:xsi=   xsi:   xsi:schemaLocation= &gt;</pre>	<pre>version="1.0" xmlns="{this-- namespace-URI- }" xmlns:xsi="{U- RL-of-governi- ng-XML-Schema- }" xsi:schemaLoc- ation="{this-- namespace-URI- }" {URL-to-find-- XML-Schema}"</pre> <pre>version="1.0" cansas1d/1.0" ="http://www.w3 schemaLocation= canSAS/1dwg/tru</pre>
---------	-----------	--------	---	--

2.2 SASroot element

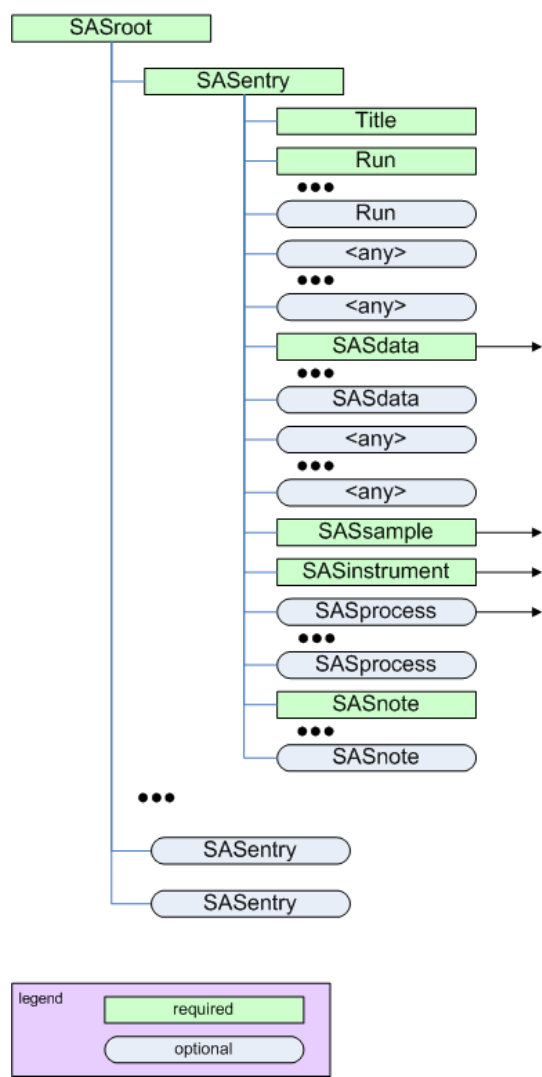

The SASroot element

Figure 2.2: The SASroot element

- parent: XML header

Table 2.2:

Name	Type	occurrence	Description	Attributes
------	------	------------	-------------	------------

Table 2.2: (continued)

<b>SASentry</b>	container	[1..unbounded]	A single SAS scan is reported in a <b>SASentry</b> . Include as many <b>SASentry</b> elements as desired. They may contain related or unrelated data. <b>name</b> is an optional attribute to provide a string for this <b>SASentry</b> . (Use of this string is not defined by this standard.)	name=short-name
-----------------	-----------	----------------	---	-----------------

2.3 SASentry element

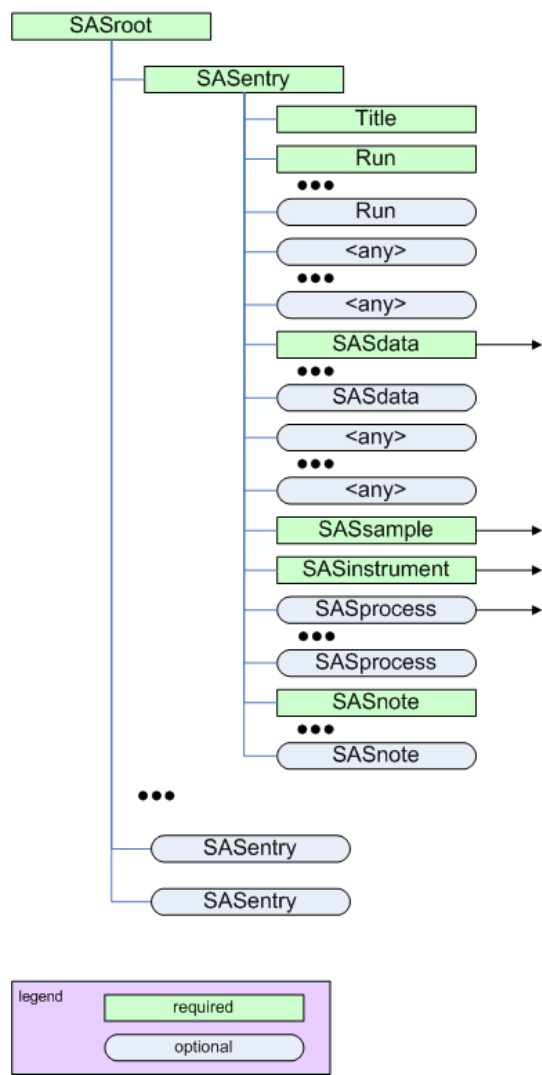


Figure 2.3: The SASentry element

- parent: SASroot

Table 2.3:

Name	Type	occurrence	Description	Attributes
<b>Title</b>	string	[1..1]	Title of this <b>SASentry</b> .	




Table 2.3: (continued)

<b>Run</b>	string	[1..unbounded]	Run identification for this <code>""SASentry""</code> . For many facilities, this is an integer. Use multiple instances of <b>Run</b> as needed. Note: How to correlate this with <code>""SASdata""</code> and <code>""SASinstrument""</code> configurations has not yet been defined. <b>name</b> is an optional string attribute to identify this particular <b>Run</b> . Could use this to associate (correlate) multiple <code>""SASdata""</code> elements with <b>Run</b> elements. (Give them the same <i>short-Run-identifier</i> .)	<b>name</b> = { short-Run-identifier }
<code>"""" { any } """"</code>	container	[0..unbounded]	Any element(s) not defined in the cansas1d/1.0 standard can be placed at this point. See <code>"" { any } ""</code> for more details.	<i>xmlns: { foreign-prefix } </i> = { foreign-namespace }
<code>""SASdata""</code>	container	[1..unbounded]	Reduced 1-D SAS data for this <code>""SASentry""</code> . Use multiple <code>""SASdata""</code> elements to represent multiple frames. Use this to associate (correlate) multiple <code>""SASdata""</code> elements with <b>Run</b> elements. (Give them the same name.)	<b>name</b> = { short-Run-identifier }
<code>"""" { any } """"</code>	container	[0..unbounded]	Any element(s) not defined in the cansas1d/1.0 standard can be placed at this point. See <code>"" { any } ""</code> for more details.	<i>xmlns: { foreign-prefix } </i> = { foreign-namespace }
<code>""SASsample""</code>	container	[1..1]	Description of the sample.	<b>name</b> = { short-SASsample-identifier }
<code>""SASinstrument""</code>	container	[1..1]	Description of the instrument	

Table 2.3: (continued)

""SASprocess""	container	[0..unbounded]	Description of a processing or analysis step.	<i>name</i> ={ short-SASprocess-identifier }
""SASnote""	container	[1..unbounded]	Free form description of anything not covered by other elements.	<i>name</i> ={ short-SASnote-identifier }

2.4 SASdata element

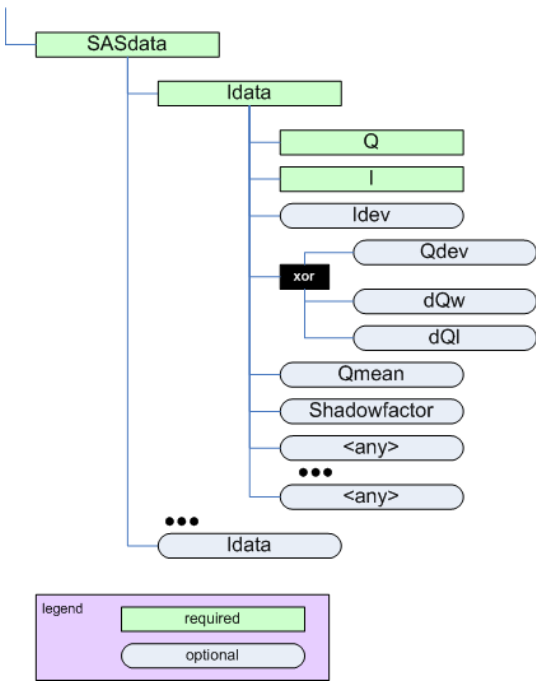

The SASdata element

Figure 2.4: The SASdata element

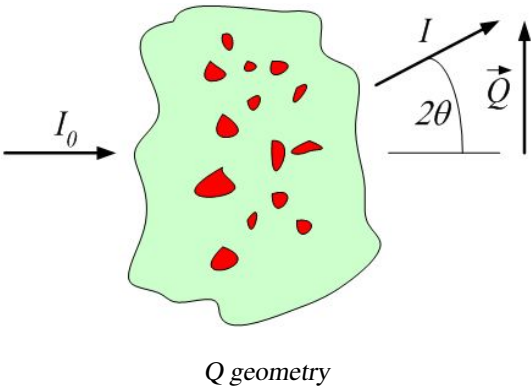


Figure 2.5: Q geometry

- parent: SASentry

2.4.1 SASdata

Table 2.4:

Name	Type	occurrence	Description	Attributes
""Idata""	container	[1..unbounded]	Idata describes a single SAS data point.	

2.4.2 Idata

Table 2.5:

Name	Type	occurrence	Description	Attributes
Q	floating-point number	[1..1]	$Q=(4\pi/\lambda)\sin(\theta)$ where $\lambda$ is the wavelength of the radiation and $2\theta$ is the angle through which the detected radiation has been scattered. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values. Either $1/\text{\AA}$ or $1/\text{nm}$ are typical.	<b>unit</b> = {units}

Table 2.5: (continued)

<b>I</b>	floating-point number	[1..1]	Intensity of the detected radiation. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values. One possibility might be <i>1/cm</i> for absolute units when the intensity describes a <i>differential cross-section per unit volume per unit solid angle</i> . <b>NOTE:</b> Be aware that there are different types of intensity that may be reported (see <a href="#">Talk:cansas1d_documentation#The_Intensity_Problem</a> ) and one should be very careful to inspect the <b>unit</b> attribute to determine how to handle subsequent data processing, especially in the area of units conversion.	<b>unit</b> = {units}
<b>Idev</b>	floating-point number	[0..1]	Estimated standard deviation of <b>I</b> . The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values. One possibility might be <i>1/cm</i> .	<b>unit</b> = {units}
<b>Qdev</b>	floating-point number	[0..1]	Estimated standard deviation of <b>Q</b> . (optional: see note below on usage)The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values. Either <i>1/A</i> or <i>1/nm</i> are typical.	<b>unit</b> = {units}

Table 2.5: (continued)

<b><i>dQ<sub>w</sub></i></b>	floating-point number	[0..1]	<b>Q</b> resolution along the axis of scanning (the high-resolution <i>slit width</i> direction). Useful for defining resolution data from slit-smearing instruments such as Bonse-Hart geometry. (optional: see note below on usage)The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values. Either <i>1/Å</i> or <i>1/nm</i> are typical.	<b>unit</b> = { units }
<b><i>dQ<sub>l</sub></i></b>	floating-point number	[0..1]	<b>Q</b> resolution perpendicular to the axis of scanning (the low-resolution <i>slit length</i> direction). Useful for defining resolution data from slit-smearing instruments such as Bonse-Hart geometry. (optional: see note below on usage)The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values. Either <i>1/Å</i> or <i>1/nm</i> are typical.	<b>unit</b> = { units }
<b><i>Q<sub>mean</sub></i></b>	floating-point number	[0..1]	Mean value of <b>Q</b> for this datum. Useful when describing data that has been binned from higher-resolution or from area detectors.The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values. Either <i>1/Å</i> or <i>1/nm</i> are typical.	<b>unit</b> = { units }

Table 2.5: (continued)

<i>Shadowfactor</i>	floating-point number	[0..1]	Describes the adjustment due to the beam stop penumbra. <b>(This definition needs revision. NIST?)</b> NOTE: There is no unit attribute.	
"" {any} ""	container	[0..unbounded]	Any element(s) not defined in the cansas1d/1.0 standard can be placed at this point. See ""{any}"" for more details.	<i>xmlns:{foreign-prefix}={foreign-namespace}</i>

### 2.4.3 Notes

1. When an optional element (*Idev*, *Qdev*, ...) is used, it must be given in every ""Idata"" within the enclosing ""SASdata"".
2. If either *dQw* or *dQl* are used, then *Qdev* is not permitted to be used.


---

#### Note

**(The Shadowfactor attribute definition needs revision. NIST?)**

---

## 2.5 SASsample element

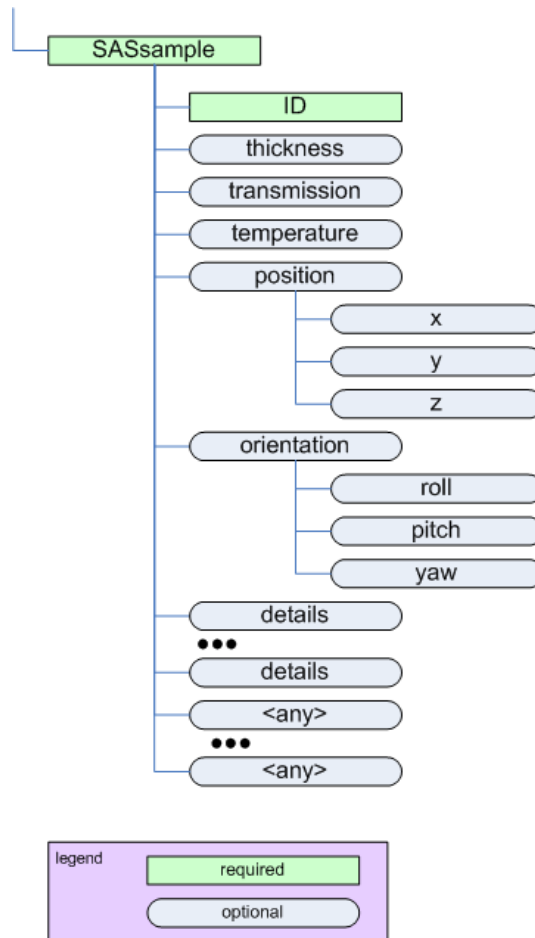

The SASsample element

Figure 2.6: The SASsample element

- parent: SASentry

Table 2.6:


Name	Type	occurrence	Description	Attributes
<b>ID</b>	string	[1..1]	Text string that identifies this sample.	
<i>thickness</i>	floating-point number	[0..1]	Thickness of this sample. Must specify the unit as an attribute.	<b>unit</b> = { unit }
<i>transmission</i>	floating-point number	[0..1]	Transmission (1-attenuation) of this sample. Express this as a fraction, not as a percentage. NOTE: there is no unit attribute.	

Table 2.6: (continued)

<i>temperature</i>	floating-point number	[0..1]	Temperature of this sample. Must specify the unit as an attribute.	<b>unit</b> = { unit }
""position""	container	[0..1]	Location in X, Y, and Z of the sample.	
""orientation""	container	[0..1]	Orientation (rotation) of the sample.	
<i>details</i>	string	[0..unbounded]	Any additional sample details.	
""{ any }""	container	[0..unbounded]	Any element(s) not defined in the cansas1d/1.0 standard can be placed at this point. See ""{ any }"" for more details.	<i>xmlns:{foreign-prefix}</i> = { foreign-namespace }

## 2.5.1 geometry

Table 2.7:

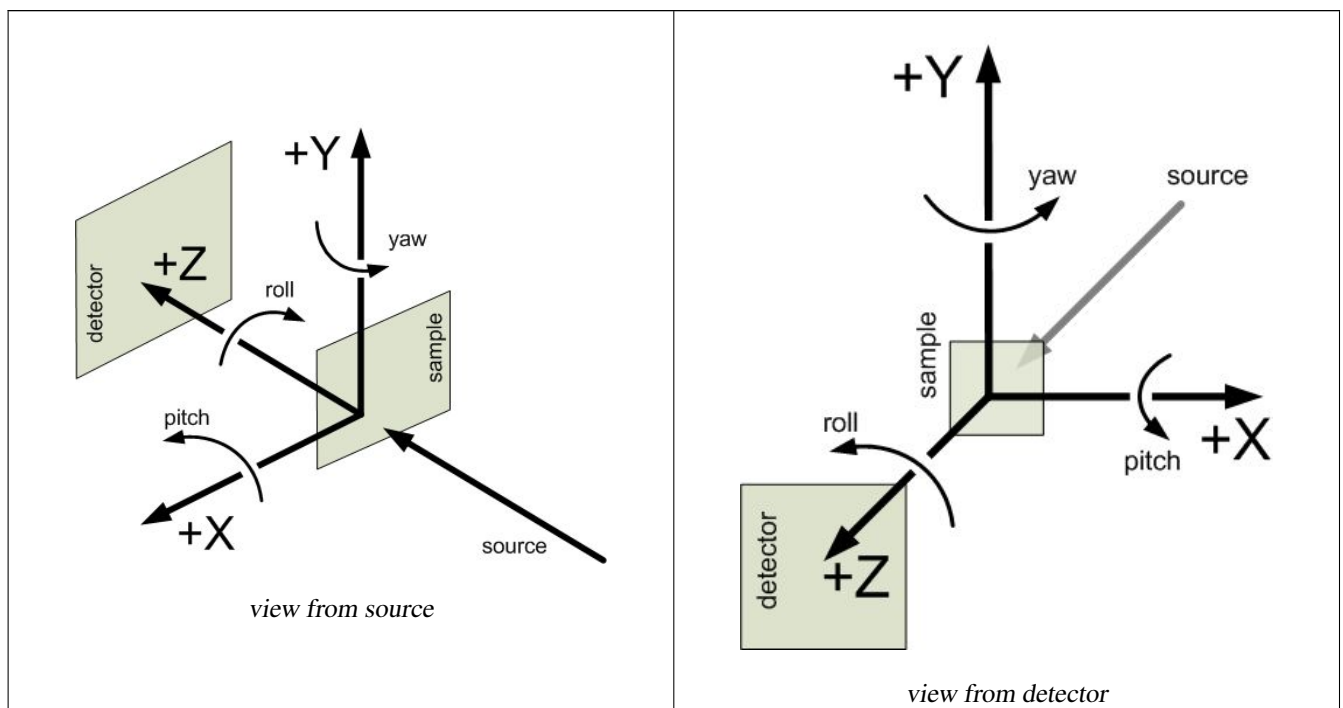

## 2.5.2 position

Table 2.8:

Name	Type	occurrence	Description	Attributes
------	------	------------	-------------	------------



Table 2.8: (continued)

<b>x</b>	floating-point number	[0..1]	Position of the sample in X. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = { units }
<b>y</b>	floating-point number	[0..1]	Position of the sample in Y. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = { units }
<b>z</b>	floating-point number	[0..1]	Position of the sample in Z. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values. Note: While Z dimension is allowed by the standard (provided by use of a standard element in the XML Schema), it does not make sense for small-angle scattering.	<b>unit</b> = { units }

### 2.5.3 orientation


Table 2.9:

Name	Type	occurrence	Description	Attributes
<b>roll</b>	floating-point number	[0..1]	Rotation about the Z axis (roll). The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = { units }
<b>pitch</b>	floating-point number	[0..1]	Rotation about the X axis (pitch). The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = { units }
<b>yaw</b>	floating-point number	[0..1]	Rotation about the Y axis (yaw). The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = { units }

**Note**

The **orientation** element is intended to describe simple rotations about a single axis rather than a full set of rotations as in a crystallographic context.

## 2.6 SASinstrument element

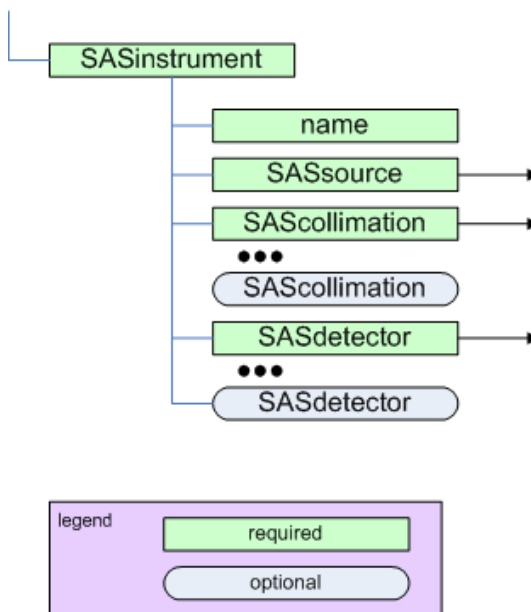

*The SASinstrument element*


Figure 2.7: The SASinstrument element

- parent: SASentry

Table 2.10:

Name	Type	occurrence	Description	Attributes
<b>name</b>	string	[1..1]	Text string that identifies the name of this instrument.	
""SASsource""	container	[1..1]	Text string that identifies the name of this source of radiation.	<b>name</b> ={name}
""SAScollimation""	container	[1..unbounded]	Text string that identifies the name of this instrument collimation.	<b>name</b> ={name}
""SASdetector""	container	[1..unbounded]	Text string that identifies the name of this detector.	

2.7 SASsource element

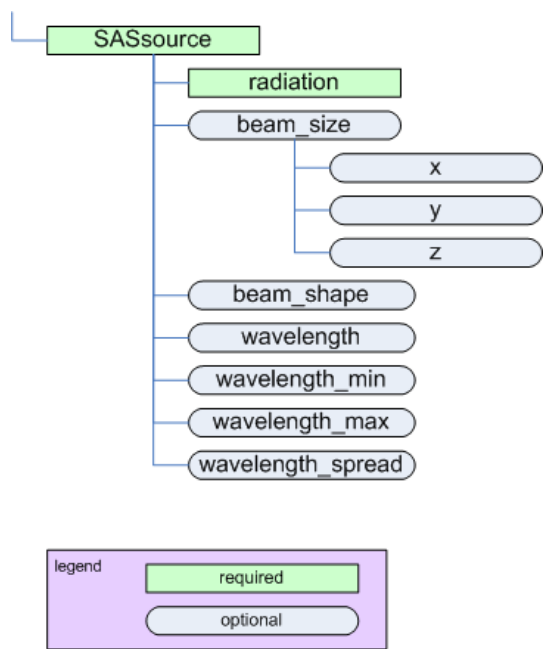

The SASsource element

Figure 2.8: The SASsource element

- parent: SASinstrument

2.7.1 SASsource

Table 2.11:

Name	Type	occurrence	Description	Attributes
------	------	------------	-------------	------------

Table 2.11: (continued)

<b>radiation</b>	string	[1..1]	Name of the radiation used. For maximum compatibility with NeXus, use one of the names defined by either <i>NeXus</i> <i>""NXsource""/type</i> * Spallation Neutron Source * Pulsed Reactor Neutron Source * Reactor Neutron Source * Synchrotron X-ray Source * Pulsed Muon Source * Rotating Anode X-ray * Fixed Tube X-ray or <i>NeXus</i> <i>""NXsource""/probe</i> * neutron * x-ray * muon * electron	
""beam_size""	container	[0..1]	Physical dimension of the beam (incident on the sample). Note: If beam is round, just use X dimension. Note: While Z dimension is allowed by the standard, it does not make sense for small-angle scattering.	<i>name</i> ={ name }
<i>beam_shape</i>	string	[0..1]	Text description of the shape of the beam (incident on the sample).	
<i>wavelength</i>	floating-point number	[0..1]	wavelength (lambda;) of radiation incident on the sample.	<b>unit</b> ={ unit }
<i>wavelength_min</i>	floating-point number	[0..1]	Some facilities specify wavelength using a range. The minimum of such a range is given by wavelength_min.	<b>unit</b> ={ unit }
<i>wavelength_max</i>	floating-point number	[0..1]	Some facilities specify wavelength using a range. The maximum of such a range is given by wavelength_max.	<b>unit</b> ={ unit }

Table 2.11: (continued)

<i>wavelength_spread</i>	floating-point number	[0..1]	Some facilities specify the width of the wavelength spectrum. The width of such a range is given by wavelength_spread.	<b>unit</b> = { unit }
--------------------------	-----------------------	--------	--	------------------------

2.7.2 beam\_size

Table 2.12:

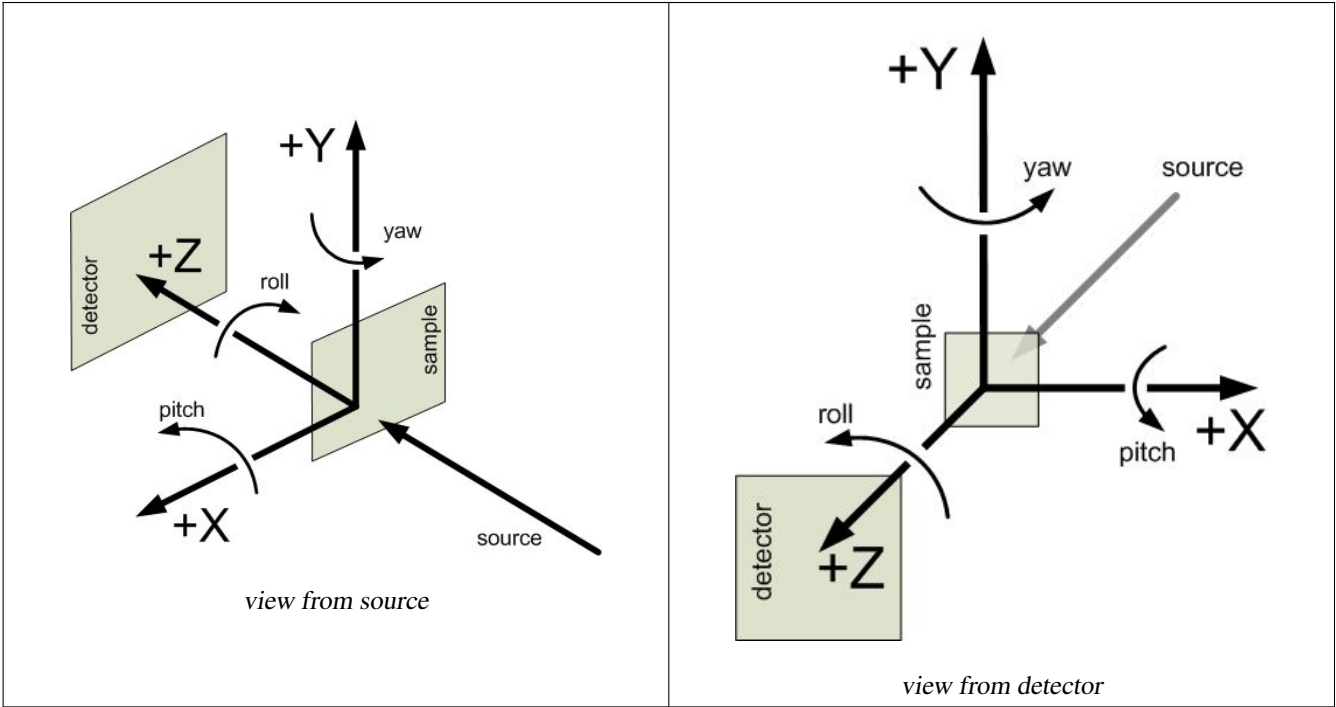



Table 2.13:

Name	Type	occurrence	Description	Attributes
<b>x</b>	floating-point number	[0..1]	Dimension of the beam in X. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = { units }
<b>y</b>	floating-point number	[0..1]	Dimension of the beam in Y. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = { units }

Table 2.13: (continued)

<b>z</b>	floating-point number	[0..1]	Dimension of the beam in Z. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values. Note: While Z dimension is allowed by the standard (provided by use of a standard element in the XML Schema), it does not make sense for small-angle scattering.	<b>unit</b> = { units }
----------	-----------------------	--------	---	-------------------------

2.8 SAScollimation element

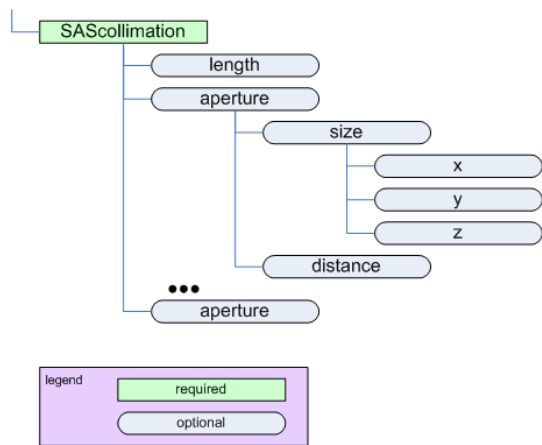

The SAScollimation element

Figure 2.9: The SAScollimation element

- parent: SASinstrument

2.8.1 SAScollimation

Table 2.14:

Name	Type	occurrence	Description	Attributes
<i>length</i>	floating-point number	[0..1]	Amount/length of collimation inserted (on a SANS instrument)	<b>unit</b> = { unit }

Table 2.14: (continued)

""aperture""	container	[0..unbounded]	Description of a slit or aperture. <b>name</b> : Optional name attribute for this aperture. <b>type</b> : Optional text attribute to describe the type aperture (pinhole, 4-blade slit, Soller slit, ...).	<b>name</b> = { name } <b>type</b> = { type }
--------------	-----------	----------------	--	--

2.8.2 aperture

Table 2.15:

Name	Type	occurrence	Description	Attributes
""size""	container	[0..1]	Opening dimensions of this aperture.	<b>name</b> = { name }
<i>distance</i>	floating-point number	[0..1]	Distance from this collimation element to the sample.	<b>unit</b> = { unit }

2.8.3 size

Table 2.16:


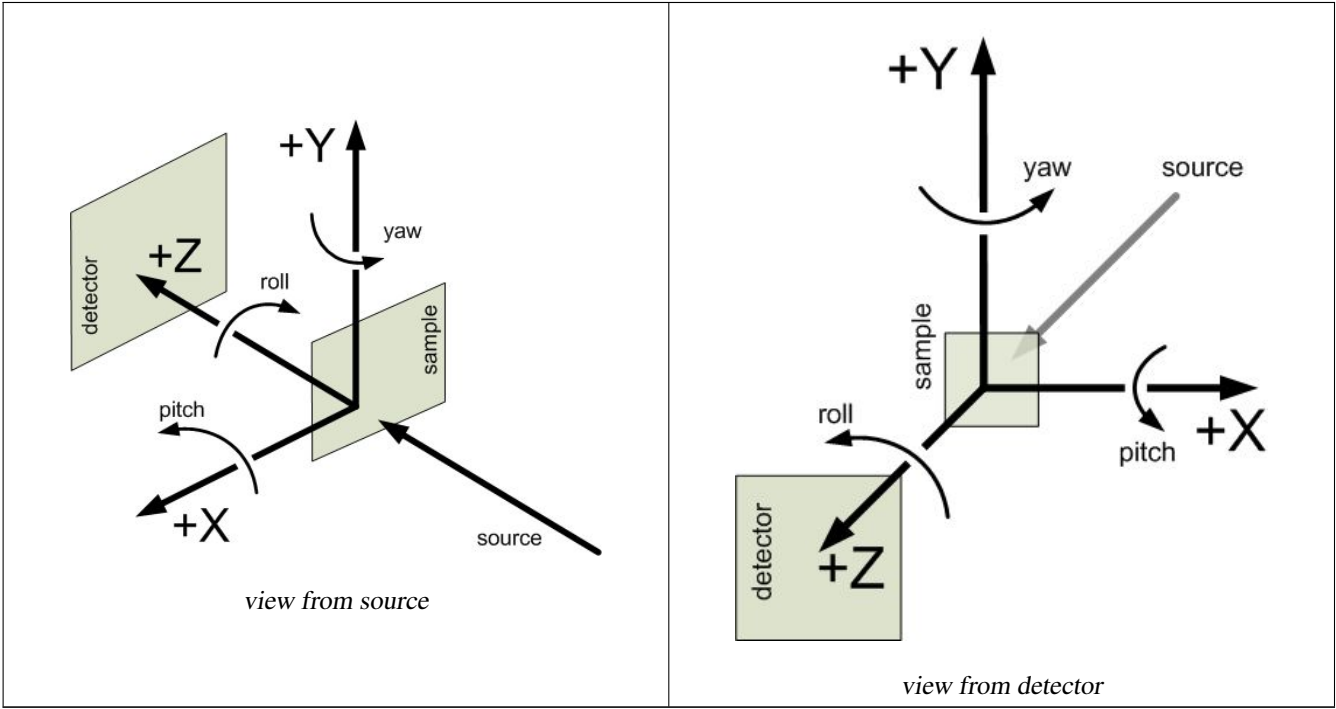


Table 2.17:

Name	Type	occurrence	Description	Attributes
<b>x</b>	floating-point number	[0..1]	Dimension of the collimation in X. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = {units}
<b>y</b>	floating-point number	[0..1]	Dimension of the collimation in Y. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = {units}
<b>z</b>	floating-point number	[0..1]	Dimension of the collimation in Z. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values. Note: While Z dimension is allowed by the standard (provided by use of a standard element in the XML Schema), it does not make sense for small-angle scattering.	<b>unit</b> = {units}



## 2.9 SASdetector element

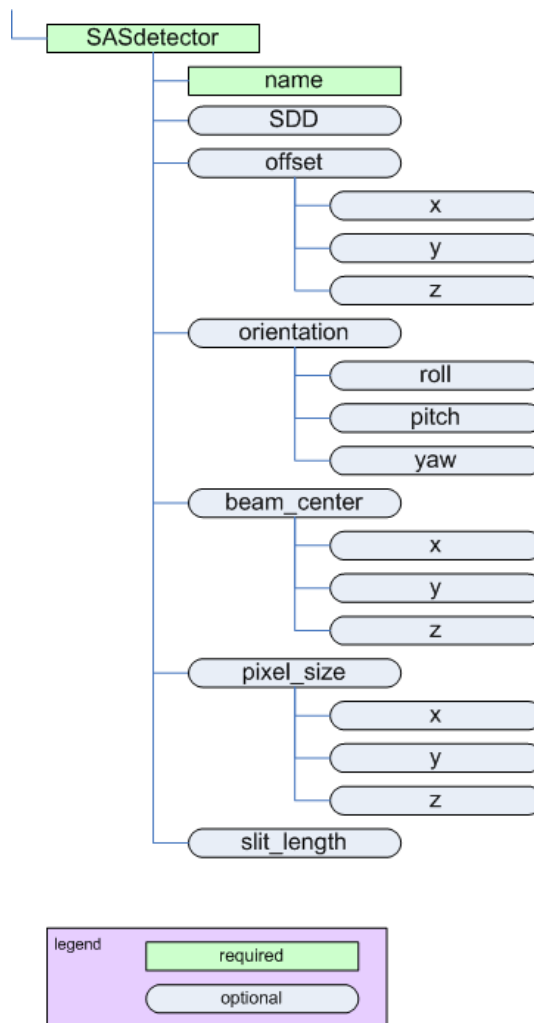

*The SASdetector element*

Figure 2.10: The SASdetector element

- parent: SASinstrument

### 2.9.1 SASdetector

Table 2.18:


Name	Type	occurrence	Description	Attributes
<b>name</b>	string	[1..1]	Text string that identifies the name of this detector.	
<b>SDD</b>	floating-point number	[0..1]	Distance between sample and detector.	<b>unit</b> = {unit}
"" offset ""	container	[0..1]	Offset of this detector position in X, Y, (and Z if necessary).	

Table 2.18: (continued)

""orientation""	container	[0..1]	Orientation (rotation) of this detector in roll, pitch, and yaw.	
""beam_center""	container	[0..1]	Center of the beam on the detector in X and Y (and Z if necessary).	
""pixel_size""	container	[0..1]	Size of detector pixels in X and Y (and Z if necessary).	
<i>slit_length</i>	floating-point number	[0..1]	Slit length of the instrument for this detector. This is expressed in the same units as Q (reciprocal space units).	<i>unit</i> = { unit }

## 2.9.2 geometry

Table 2.19:

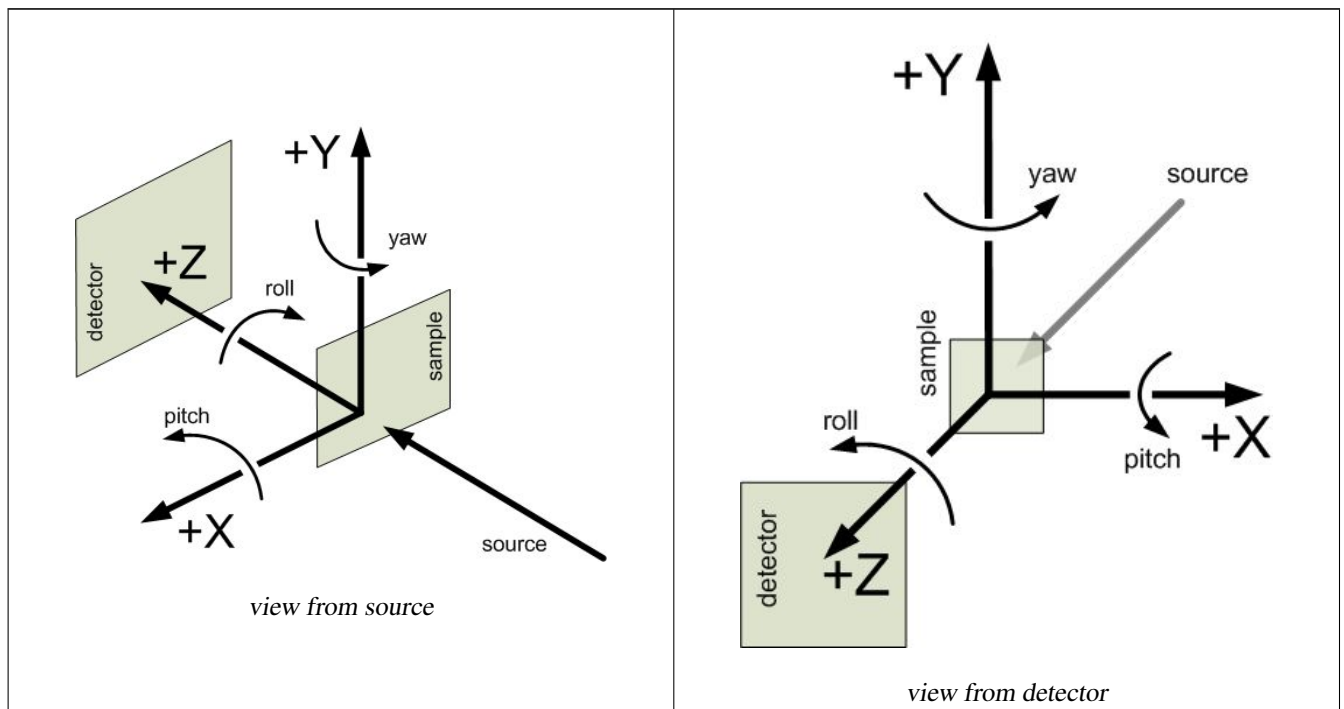

## 2.9.3 offset

Table 2.20:

Name	Type	occurrence	Description	Attributes
------	------	------------	-------------	------------

Table 2.20: (continued)

<b>x</b>	floating-point number	[0..1]	Offset of the detector position in X. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = { units }
<b>y</b>	floating-point number	[0..1]	Offset of the detector position in Y. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = { units }
<b>z</b>	floating-point number	[0..1]	Offset of the detector position in Z. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values. Note: While Z is allowed by the standard (provided by use of a standard element in the XML Schema), it does not make sense for small-angle scattering.	<b>unit</b> = { units }

## 2.9.4 orientation

Table 2.21:

Name	Type	occurrence	Description	Attributes
<b>roll</b>	floating-point number	[0..1]	Rotation about the Z axis (roll). The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = { units }
<b>pitch</b>	floating-point number	[0..1]	Rotation about the X axis (pitch). The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = { units }
<b>yaw</b>	floating-point number	[0..1]	Rotation about the Y axis (yaw). The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = { units }

## 2.9.5 beam\_center

Table 2.22:

Name	Type	occurrence	Description	Attributes
<b>x</b>	floating-point number	[0..1]	Position of the beam center on the detector in X. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = {units}
<b>y</b>	floating-point number	[0..1]	Position of the beam center on the detector in Y. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = {units}
<b>z</b>	floating-point number	[0..1]	Position of the beam center on the detector in Z. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values. Note: While Z dimension is allowed by the standard (provided by use of a standard element in the XML Schema), it does not make sense for small-angle scattering.	<b>unit</b> = {units}

## 2.9.6 pixel\_size


Table 2.23:

Name	Type	occurrence	Description	Attributes
<b>x</b>	floating-point number	[0..1]	Size of a detector pixel in X. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = {units}

Table 2.23: (continued)

y	floating-point number	[0..1]	Size of a detector pixel in Y. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values.	<b>unit</b> = { units }
z	floating-point number	[0..1]	Size of a detector pixel in Z. The <b>unit</b> attribute is required. See <a href="#">cansas1d_documentation#Rules</a> for acceptable values. Note: While Z dimension is allowed by the standard (provided by use of a standard element in the XML Schema), it does not make sense for small-angle scattering.	<b>unit</b> = { units }

2.10 SASprocess element

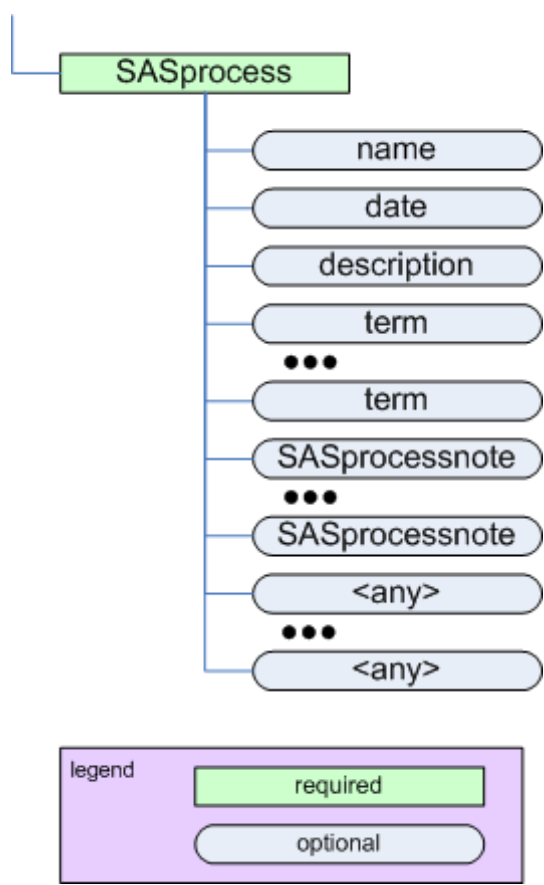

The SASprocess element

Figure 2.11: The SASprocess element

- parent: SASentry

Table 2.24:

Name	Type	occurrence	Description	Attributes
<i>name</i>	string	[0..1]	Optional name for this data processing or analysis step.	
<i>date</i>	string	[0..1]	Optional date for this data processing or analysis step. Use a format which is easily machine-readable such as yyyy-mm-dd hh:mm:ss The format for the date string may be specified at a later date.	

Table 2.24: (continued)


<i>description</i>	string	[0..1]	Optional description for this data processing or analysis step.	
<i>term</i>	string	[0..unbounded]	This is used to specify the value of a single variable, parameter, or term (while defined here as a string, it could be a number) related to the SASprocess step.	<b>unit</b> = { unit }
""SASprocessnote""	container	[1..unbounded]	This element is used to describe anything about SASprocess that is not already described.	
"" { any } ""	container	[0..unbounded]	Any element(s) not defined in the cansas1d/1.0 standard can be placed at this point. See "" { any } "" for more details.	<i>xmlns: { foreign-prefix } </i> = { foreign-namespace }

### 2.10.1 SASprocessnote

Table 2.25:

Name	Type	occurrence	Description	Attributes
"" { any } ""	container	[0..unbounded]	Any element(s) not defined in the cansas1d/1.0 standard can be placed at this point. See "" { any } "" for more details.	<i>xmlns: { foreign-prefix } </i> = { foreign-namespace }

2.11 SASnote element

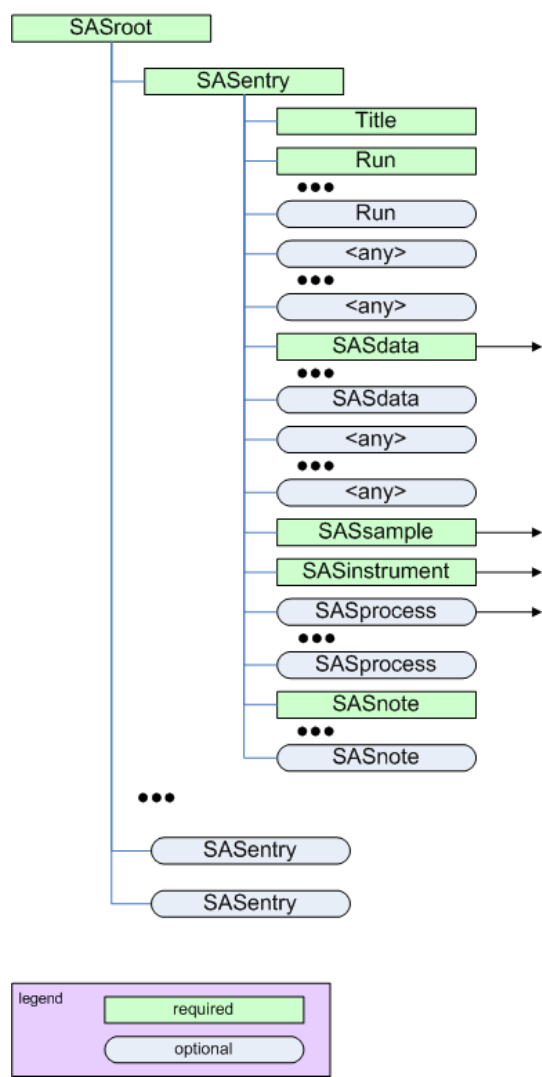

The SASroot element

Figure 2.12: The SASroot element

- parent: SASentry

Table 2.26:

Name	Type	occurrence	Description	Attributes
"" {any} ""	container	[0..unbounded]	Any element(s) not defined in the cansas1d/1.0 standard can be placed at this point. See ""{any}"" for more details.	<i>xmlns:{foreign-prefix}</i> ={foreign-namespace}



## 2.12 "{any}" element

Table 2.27:

Name	Type	occurrence	Description	Attributes
<b>{any}</b>	container	[0..unbounded]	Any element(s) not defined in the cansas1d/1.0 standard can be placed at this point. (These are called <i>foreign</i> elements. It is suggested to associate foreign elements with a foreign namespace to differentiate them from the canSAS elements in the XML file.)	<b><i>xmlns:{foreign-prefix}</i></b> ={ foreign-namespace }

## Chapter 3

# Bindings and Software Support

Bindings (import/export drivers) and other software support have been created and contributed. These are listed here by the language or software environment.

### 3.1 Fortran binding

The development of the FORTRAN language, so beloved of scientists, pre-dates the development of XML. And it shows. FORTRAN is not a language that manipulates strings with ease, and this makes parsing XML decidedly awkward. So unless you *really* have to use FORTRAN, you are probably better off with C/C++ (or something else more 'modern'), see for example Daniel Veillard's LIBXML2 library at <http://xmlsoft.org/> or Frank van den Berghen's parser at <http://www.applied-mathematics.net/tools/xmlParser.html>.

If you have to use a dialect earlier than FORTRAN-90 (F90), then the chances are you will have to code your own parser.

#### 3.1.1 Software Development Kits

For later dialects, there are some SDK's available on the Web:

- F90:
  - XMLPARSE - by Arjen Markus at <http://xml-fortran.sourceforge.net/>
  - FoX - by Toby White others at <http://uszla.me.uk/space/software/FoX/>
- For F95:
  - XML - by Mart Rentmeester at <http://nn-online.org/code/xml/>

#### 3.1.2 canSAS 1-D SAS XML v1.0 support

Steve King[mailto:s.m.king@rl.ac.uk] (ISIS) has provided a F77 routine (*SASXML\_G77.F*) that will read CanSAS XML v1.0 files.

### 3.2 IgorPro binding

An import tool (binding) for IgorPro has been created (cansasXML.ipf). You can check out the IgorPro working directory from the SVN server (see instructions below).

To use the canSASxml.ipf procedure, you must have the XMLutils XOP IGOR plugin installed. See the Usage Notes below.

---

---

**Note**

Note that this tool is not a true **binding** in that the structure of the XML file is not replicated in IgorPro data structures. This tool reads the vectors of 1-D SAS data (Q, I, ...) into IgorPro waves (Qsas, Isas, ...). The tool also reads most of the metadata into an IgorPro textWave for use by other support in IgorPro.

---

**Note**

Note that the code described here is *not a complete user interface*. (See further comments below.) It is expected that this code will be called by a graphical user interface routine and that routine will handle the work of copying the loaded SAS data in IgorPro from the root:Packages:CS\_XMLreader data folder to the destination of choice (including any renaming of waves as desired).

---

**file**

*cansasXML.ipf*

**author**

Pete R. Jemian <jemian@anl.gov>

**date**

2009-09-02

**version**

1.11 (**requires** latest XMLutils XOP -- see below)

**purpose**

implement an IgorPro file reader to read the canSAS 1-D reduced SAS data in XML files adheres to the cansas1d/1.0 standard: [http://www.smallangles.net/wgwiki/index.php/cansas1d\\_documentation](http://www.smallangles.net/wgwiki/index.php/cansas1d_documentation)

**URL**

TRAC: <http://svn.smallangles.net/trac/canSAS/browser/1dwg/trunk/IgorPro/cansasXML.ipf>

SVN: <http://svn.smallangles.net/svn/canSAS/1dwg/trunk/IgorPro/cansasXML.ipf>

**requires**

IgorPro: <http://www.wavemetrics.com> XMLutils - XOP: <http://www.igorexchange.com/project/XMLutils> (IGOR.5.04.x-1.x-dev, 2008-Aug-22)

### 3.2.1 Checkout of support code in Subversion

Subversion (<http://subversion.tigris.org/>) is a program for managing software versions. There are command line and GUI clients for a variety of operating systems. We won't recommend any here but will show the command lines necessary.

#### 3.2.1.1 XMLutils XOP

The XMLutils XOP, written by Andrew Nelson (ANSTO), is hosted on the IgorExchange (<http://www.igorexchange.com/>).

One good location to place the checked out XMLutils directory is in the Wavemetrics directory, next to the Igor Pro Folder

```
svn://svn.igorexchange.com/packages/XMLutils/ XMLutils
```

In the future, to retrieve an updated version of this support, go into the XMLutils directory (created above) and type the command

```
svn update
```

This will check the repository and update files as needed. If the installer program was updated, you'll need to run the new installer program. It is not necessary to uninstall first.

The installer executables contained in the download will do all the installation for you. They will place the XOP in the folder */User Procedures/motofit/XMLutils*, and create a shortcut/alias to the plugin in */Igor Extensions*. Packages from other facilities should place the XOP there as well.

---

### 3.2.1.2 cansasXML.ipf

Check out the canSAS 1d SAS XML reader from the subversion repository:

```
svn checkout http://svn.smallangles.net/svn/canSAS/ldwg/trunk cansas-ldwg
```

This will download lots of extra files. The file of interest is in the IgorProdirectory and is called cansasXML.ipf

In the future, to retrieve an updated version of this support, go intothe cansas-ldwg directory (created above) and type the command

```
svn update
```

This will check the repository and update files as needed.

## 3.2.2 Installation

1. License and Install IgorPro (should have already been done by now)
2. Quit IgorPro if it is running
3. Download XMLutils XOP. Either checkout from subversion (see above) or, with a web browser, go to <http://svn.igorexchange.com/>
4. Install XMLutils XOP by double-clicking the installer for you operating system.
5. Download cansasXML.ipf. Either checkout from subversion (see above) or, with a web browser, copy cansasXML.ipf from on-line subversion repository: <http://svn.smallangles.net/svn/canSAS/ldwg/trunk/IgorPro/cansasXML.ipf>
6. Copy cansasXML.ipf file to ...WavemetricsIgor Pro FolderUser Procedures (or file system equivalent)
7. Then, you should be able to restart IgorPro and progress from there

## 3.2.3 Usage Notes

To use the canSASxml.ipf procedure, you must have the XMLutils XOP IGOR plugin installed. This may be downloaded from the IgorExchange Project site. There are installer executables contained in the download that will do all the installation for you. Each installer will place the XOP in the folder ...WavemetricsIgor Pro FolderUser ProceduresmotofitXMLutils, and create a shortcut/alias to the plugin in ...WavemetricsIgor Pro FolderIgor Extensions.

## 3.2.4 What it does

Given an XML file, **CS\_XmlReader(fileName)** attempts to open the fileand read its contents as if it conformed to the canSAS XML standard forreduced 1-D SAS data (cansas1d/1.0, also known as SASXML). If the file is found to be non-conforming, then **CS\_XmlReader(fileName)** returns with anerror code (show below), otherwise it returns **0** that indicates *no error*.All data read by this code is left in the IgorPro data folder **root:Packages:CS\_XMLreader** for pickup by the calling routine.(Two examples are provided to show how a routine might retrieve the data.)

After opening the XML file (with a file identifier *fileID*), control is passed to **CS\_1i\_parseXml(fileID)** which then walks through the XMLelements. For each **SASentry** in the file, a new data folder is createdwith the name derived from the Title element (or best effort determination).Efforts are taken to avoid duplication of data folder names (using standard IgorPro routines). For **SASentry** elements that contain more than one SASdataelement, a **SASdata** folder is created for each and the corresponding *I(Q)* is placed in that subfolder. When only one **SASdata** is found, the *I(Q)*data is placed in the main *Title* folder.

### 3.2.4.1 data columns

Each column of data in the **SASdata/Idata/\*** table is placed into a singleIgorPro wave. At present, the code does not check for non-standard data columns.(The capability is built into the code but is deactivated at present).

### 3.2.4.2 metadata

Additional metadata is collected into a single text wave (*metadata*) where the first column is an identifier (or *key*) and the second identifier is the *value*. Only those keys with non-empty values are retained in the metadata table. **CAUTION:** The *values* are not checked for characters that may cause trouble when placed in a wave note. This will be the responsibility of the calling routine to *clean these up* if the need arises.

The code checks for most metadata elements and will check for repeated elements where the standard permits.

Here is an example of the metadata for the "cs\_collagen\_full.xml" case study:

Table 3.1: metadata for the "cs\_collagen\_full.xml" case study

<i>i</i> (row)	metadata[ <i>i</i> ][0] ( <i>key</i> )	metadata[ <i>i</i> ][1] ( <i>value</i> )
0	xmlFile	cs_collagen_full.xml
1	namespace	cansas1d/1.0
2	Title	dry chick collagen, d = 673 A, 6531 eV, X6B
3	Run	Sep 19 1994 01:41:02 am
4	SASsample/ID	dry chick collagen, d = 673 A, 6531 eV, X6B
5	SASinstrument/name	X6B, NSLS, BNL
6	SASinstrument/SASsource/radiation	X-ray synchrotron
7	SASinstrument/SASsource/wavelength	1.898
8	SASinstrument/SASsource/wavelength/@Anit	
9	SASinstrument/SASdetector/@name	X6B PSD
10	SASnote	Sep 19 1994      01:41:02   ↔ am      Elt: 00090 Seconds ID: No spectrum   ↔ identifier defined Memory Size: 8192 Chls   ↔ Conversion Gain: 1024   Adc Offset  dry chick collagen, d =   ↔ 673 A 6531 eV, X6B

### 3.2.4.3 XML foreign namespace elements

These are ignored at this time.

### 3.2.4.4 XML namespace and header

The routine does a *best-efforts* check to ensure that the given XML file conforms to the [[cansas1d\_documentation#Required\_XML\_file\_header|required XML file header]]. If you take a minimalist view (*a.k.a.* a shortcut), it is likely that your file may be refused by this and other readers. Pay particular attention to UPPER/lower case in the text **cansas1d/1.0** as this is a **key component** used to index through the XML file.

### 3.2.4.5 XML stylesheet processing-instruction is not generated

The "*XMLutils*" package does not provide a method to insert the prescribed XML stylesheet processing-instruction into the XML data file.

```
?xml-stylesheet type=text/xsl href=example.xsl ?
```

If this processing-instruction is desired, it must be added to each XML data file by other methods such as use of a text editor or application of an XSLT transformation.

### 3.2.5 List of Functions

These are (most of) the FUNCTIONS in the cansasXML.ipf code. The only functions of interest are **CS\_XmlReader(fileName)** which reads the named XML file and loads SAS data and the two demonstration functions **prj\_grabMyXmlData()** and **prjTest\_cansas1d()** together show a usage example.

- **CS\_XmlReader(fileName)** : open a canSAS 1-D reduced SAS XML data file
- input: *fileName* (string) name of canSAS XML file (can include file system path name to file)
- returns:
  - 0 successful
  - -1: XML file not found
  - -2: root element is not SASroot with valid canSAS namespace
  - -3: SASroot version is not 1.0
  - -4: no SASentry elements (NOT USED NOW)
  - -5: XOPutils needs upgrade
- **CS\_li\_parseXml(fileID)** : **This is what guides the work**, given a file ID returned from **XMLOpenFile()**, parses that file for SAS data and metadata (li in the function name signifies this is a function that supports INPUT from version 1.0 XML files)
- **CS\_li\_getOneSASdata(fileID, Title, SASdataPath)** : harvest the data and metadata in the specific SASdata element
- **CS\_li\_getOneVector(file, prefix, XML\_name, Igor\_name)** : harvest just one column (vector) of data
- **CS\_li\_GetReducedSASdata(fileID, SASdataPath)** : grab the data and put it in the working data folder
- **CS\_li\_locateTitle(fileID, SASentryPath)** : determine the title for this experiment
- **CS\_appendMetaData(fileID, key, xpath, value)** : queries XML file for **xpath**. If **value** is not empty, appends it to **metadata** where *last* is the new last row: `metadata[last][0]=key; metadata[last][1]=value`
- **CS\_buildXPathStr(prefix, value)** : this function can be used only with very simple XPath constructions
- **CS\_cleanFolderName(proposal)** : given a proposal string, returns a candidate folder name for immediate use
- **CS\_findElementIndex(matchStr)** : looks for element index in structure *W\_ElementList* returned from call to **XmlElementList(fileID)**
- **CS\_getDefaultNamespace(fileID)** : returns the string containing the default namespace for the XML file
- **CS\_registerNameSpaces()** : Builds a table of all namespaces used in the XML file and appends **W\_ElementList** with full namespace-xpath string for each element.
- **CS\_simpleXmlListXPath(fileID, prefix, value)** : Calls **XMLlistXPath()** with proper namespace prefix attached.
- **CS\_simpleXmlWaveFmXPath(fileID, prefix, value)** : Calls **XMLwaveFmXPath()** with proper namespace prefix attached.
- **CS\_updateWaveNote(wavName, key, value)** : adds (or replaces) definition of *key=value* in the wave note of *wavName*
- **CS\_XmlStrFmXPath(fileID, prefix, value)** : Calls **XmlStrFmXPath()** with proper namespace prefix attached. Trims the result string.
- **CS\_XPath\_NS(simpleStr)** : this function adds namespace info as necessary to simpleStr (an XPath)
- **TrimWS(str)** : Calls **TrimWSL(TrimWSR(str))**

- TrimWSL(str) : Trims white space from left (leading) end of **str**
- TrimWSR(str) : Trims white space from right (trailing) end of **str**
- prjTest\_cansas1d() : Demonstration function that calls **CS\_XmlReader(fileName)** for many of the test data sets.
- prj\_grabMyXmlData() : Demonstration function that moves loaded data from root:Packages:CS\_XMLreader to a user's data folder. (In this *example*, that folder is root:PRJ\_canSAS.)
- testCollette() : Demonstration function that reads an ISIS/LOQ file and copies the data to the root folder a la COLLETE

### 3.2.6 Example test case

Here is an example running the test routine **prjTest\_cansas1d()**.

```
•prjTest_cansas1d()
XMLopenfile: File(path) to open doesn't exist, or file can't be opened
elmo.xml either not found or cannot be opened for reading
  Completed in 0.00669666 seconds
XMLopenfile: XML file was not parseable
cansasXML.ipf: failed to parse XML
  Completed in 0.0133704 seconds
root element is not SASroot with valid canSAS namespace
  Completed in 0.0134224 seconds
bimodal-test1.xml      identified as: cansas1d/1.0 XML file
  Title: SAS bimodal test1
  Completed in 0.068654 seconds
root element is not SASroot with valid canSAS namespace
  Completed in 0.0172572 seconds
root element is not SASroot with valid canSAS namespace
  Completed in 0.0123102 seconds
root element is not SASroot with valid canSAS namespace
  Completed in 0.00930118 seconds
ISIS_SANS_Example.xml  identified as: cansas1d/1.0 XML file
  Title: standard can 12mm SANS
  Completed in 0.0410387 seconds
W1W2.xml              identified as: cansas1d/1.0 XML file
  Title: standard can 12mm SANS
  Title: TK49 standard 12mm SANS
  Completed in 0.0669074 seconds
ill_sasxml_example.xml identified as: cansas1d/1.0 XML file
  Title: ILL-D22 example: 7D1 2mm
  Completed in 0.0332752 seconds
isis_sasxml_example.xml identified as: cansas1d/1.0 XML file
  Title: LOQ TK49 Standard 12mm C9
  Completed in 0.0388868 seconds
r586.xml              identified as: cansas1d/1.0 XML file
  Title: ILL-D11 example1: 2A 5mM 0%D2O
  Completed in 0.0213737 seconds
r597.xml              identified as: cansas1d/1.0 XML file
  Title: ILL-D11 example2: 2A 5mM 0%D2O
  Completed in 0.0221894 seconds
xg009036_001.xml      identified as: cansas1d/1.0 XML file
  Title: det corrn 5m
  Completed in 0.0286721 seconds
cs_collagen.xml        identified as: cansas1d/1.0 XML file
  Title: dry chick collagen, d = 673 A, 6531 eV, X6B
  Completed in 0.0296247 seconds
cs_collagen_full.xml   identified as: cansas1d/1.0 XML file
  Title: dry chick collagen, d = 673 A, 6531 eV, X6B
  Completed in 0.0751836 seconds
cs_af1410.xml          identified as: cansas1d/1.0 XML file
```

```
Title: AF1410-10 (AF1410 steel aged 10 h)
Title: AF1410-8h (AF1410 steel aged 8 h)
Title: AF1410-qu (AF1410 steel aged 0.25 h)
Title: AF1410-cc (AF1410 steel aged 100 h)
Title: AF1410-2h (AF1410 steel aged 2 h)
Title: AF1410-50 (AF1410 steel aged 50 h)
Title: AF1410-20 (AF1410 steel aged 20 h)
Title: AF1410-5h (AF1410 steel aged 5 h)
Title: AF1410-1h (AF1410 steel aged 1 h)
Title: AF1410-hf (AF1410 steel aged 0.5 h)
Completed in 0.338425 seconds
XMLopenfile: File(path) to open doesn't exist, or file can't be opened
cansas1d-template.xml either not found or cannot be opened for reading
Completed in 0.00892823 seconds
1998spheres.xml identified as: cansas1d/1.0 XML file
Title: 255 nm PS spheres
Title: 460 nm PS spheres
Completed in 2.87649 seconds
XMLopenfile: File(path) to open doesn't exist, or file can't be opened
does-not-exist-file.xml either not found or cannot be opened for reading
Completed in 0.00404549 seconds
cs_rr_polymers.xml identified as: cansas1d/1.0 XML file
Title: Round Robin Polymer A
Title: Round Robin Polymer B
Title: Round Robin Polymer C
Title: Round Robin Polymer D
Completed in 0.0943477 seconds
s81-polyurea.xml identified as: cansas1d/1.0 XML file
Title: S7 Neat Polyurea
Completed in 0.0361616 seconds
```

### 3.2.7 Graphical User Interface

At least two groups are working on graphical user interfaces that use the canSAS 1-D SAS XML format binding to IgorPro. The GUIs are intended to be used with their suites of SAS analysis tools (and hide the details of using this support code from the user).

NOTE: There is no support yet for writing the data back into the canSAS format. Several details need to be described, and these are being collected on the discussion page for the XML format

#### 3.2.7.1 Irena tool suite

Jan Ilavsky's **Irena** tool suite (<http://usaxs.xor.aps.anl.gov/staff/ilavsky/irena.html>) for IgorPro has a GUI to load the data found in the XML file. Refer to the WWW site for more details.

## 3.3 Java JAXB binding

Documentation for the JAXB binding is spotty at this time. You can check it out with subversion:

```
svn checkout http://svn.smallangles.net/svn/canSAS/ldwg/trunk cansas-1d-standard
```

(where cansas-1d-standard is a local directory name).

### 3.3.1 JAXB

- Question : What is JAXB?



- Answer : Java Architecture for XML Binding (JAXB): <http://java.sun.com/developer/technicalArticles/Webservices/jaxb/>
- Wow! : Is it available for other languages?
- Answer : Ask Google. JAXB is for Java. One good advisory page: <http://www.devx.com/ibm/Article/20261>
- Question : How do I pull out the  $I(Q)$  data?
- Answer : see fragment below (gets data for desmearing)

### 3.3.2 JAXB\_cansas1d\_reader.java: example usage in JAVA

Here is a Java class that shows how to use the JAXB binding. Use this with any of the test data supplied with the cansas-1d-standard directory (above). By default, it shows what is in the *1998spheres.xml* example file: two samples. (You'll have to get the directory paths right until this documentation improves.)

```
/**
 *
 */
package jlake;

import java.io.File;
import java.util.List;

import javax.xml.bind.JAXBContext;
import javax.xml.bind.JAXBElement;
import javax.xml.bind.JAXBException;
import javax.xml.bind.Unmarshaller;
import cansas1d.SASdataType;
import cansas1d.SASentryType;
import cansas1d.SASrootType;
import cansas1d.SASentryType.Run;

/**
 * @author Pete Jemian
 *
 */
public class JAXB_cansas1d_reader {

    /**
     * @param args
     */
    @SuppressWarnings("unchecked")
    public static void main(String[] args) {
        JAXBContext jc;
        String xmlFile;
        // xmlFile = "cs_af1410.xml";
        xmlFile = "1998spheres.xml";
        try {
            // use the cansas1d/1.0 schema that is bound to a Java structure
            jc = JAXBContext.newInstance("cansas1d");
            Unmarshaller unmarshaller = jc.createUnmarshaller();
            // open the XML into a Java data structure
            JAXBElement<SASrootType> xmlJavaData = (JAXBElement<SASrootType>) unmarshaller
                .unmarshal(new File(xmlFile));
            // canSAS XML file is now loaded in memory
            SASrootType sasRootType = xmlJavaData.getValue();
            int numEntries = sasRootType.getSASentry().size();
            System.out.println("SASentry elements: " + numEntries);
            for( int i = 0; i < numEntries; i++ ) {
                System.out.println("SASentry");
                SASentryType entry = sasRootType.getSASentry().get(i);
```

```
System.out.printf("Title:\t%s\n", entry.getTitle());
List<SASentryType.Run> runs = entry.getRun();
System.out.printf("#Runs:\t%d\n", runs.size());
for( int j = 0; j < runs.size(); j++ ) {
    Run run = (Run) runs.get(j);
    System.out.printf("Run@name:\t%s\n", run.getName());
    System.out.printf("Run:\t%s\n", run.getValue());
}
List<SASdataType> datasets = entry.getSASdata();
System.out.printf("#SASdata:\t%d\n", entry.getSASdata().size());
for( int j = 0; j < runs.size(); j++ ) {
    SASdataType sdt = (SASdataType) datasets.get(j);
    System.out.printf("SASdata@name:\t%s\n", sdt.getName());
    System.out.printf("#points:\t%d\n", sdt.getIdata().size());
}
System.out.println();
}

System.out.println("the end.");

} catch (JAXBException e) {
    // TODO Auto-generated catch block
    e.printStackTrace();
    System.out.printf("Could not open (unmarshall) XML file: %s\n", xmlFile);
}
}
```

### 3.3.3 example: how to retrieve I(Q)

Look for the line that has **Qsas[i] = sdt.getIdata().get(i).getQ().getValue();** to see the operations to unwind the data into usable double[] vectors. Pretty straightforward.

- sdt : SASdataType object
- getIdata() : /SASdata/Idata
- get(i) : /SASdata/Idata[i]
- getQ() : /SASdata/Idata/Q
- getValue() : /SASdata/Idata/Q (value, not the unit)

#### 3.3.3.1 GetSASdata.java

```
/**
 *
 */
package jlake;

import java.io.File;
import java.util.List;

import javax.xml.bind.JAXBContext;
import javax.xml.bind.JAXBElement;
import javax.xml.bind.JAXBException;
import javax.xml.bind.Unmarshaller;
```

```
import cansas1d.SASdataType;
import cansas1d.SASdetectorType;
import cansas1d.SASentryType;
import cansas1d.SASinstrumentType;
import cansas1d.SASrootType;
import cansas1d.SASentryType.Run;

/**
 * @author Pete Jemian
 *
 */
public class GetSASdata {

    private static SASrootType sasRoot;    // SAS data (from cansas1d/1.0 XML file)
    private static double[] Qsas;         // input Q
    private static double[] Isas;         // input I (slit-smeared)
    private static double[] Idev;         // input Idev (slit-smeared)
    private static double[] Ismr;         // calculated I slit-smeared
    private static double[] Idsm;         // calculated I desmeared
    private static double[] IdsmDev;      // calculated Idev desmeared
    private static double slit_length;

    /**
     * @param xmlPropertyFile
     */
    public GetSASdata(String xmlDataFile)
    {
        // load SAS data into memory
        try {
            sasRoot = (SASrootType) loadXML("cansas1d", xmlDataFile);
        } catch (JAXBException e) {
            e.printStackTrace();
            System.out.println("ERROR: Cannot find or interpret SAS XML data file:\t" + ↵
                               xmlDataFile);
            return;
        }

        // SAS data are loaded
        // grab the SAS data to be desmeared
        int entryIndex = 0; // /SASentry[1] : unit base in XML, 0 base in Java
        int dataIndex = 0; // SASdata[1]
        int detectorIndex = 0; // SASdetector[1]
        SASentryType entry = (SASentryType) sasRoot.getSASentry().get(entryIndex);
        SASdataType sdt = (SASdataType) entry.getSASdata().get(dataIndex);
        if (sdt.getName().trim().compareTo("slit-smeared") != 0) {
            System.out.println("selected SASdata element must start: <SASdata name=\"slit-smeared ↵
                               \">");
            // throw something (an exception) here?
            return;
        }
        int numPoints = sdt.getIdata().size();
        Qsas = new double[numPoints]; // input Q
        Isas = new double[numPoints]; // input I (slit-smeared)
        Idev = new double[numPoints]; // input Idev (slit-smeared)
        for (int i = 0; i < numPoints; i++) {
            Qsas[i] = sdt.getIdata().get(i).getQ().getValue();
            Isas[i] = sdt.getIdata().get(i).getI().getValue();
            Idev[i] = sdt.getIdata().get(i).getIdev().getValue();
        }
    }
}
```

```
    Ismr = new double[numPoints];    // calculated I slit-smeared
    Idsm = new double[numPoints];    // calculated I desmeared
    IdsmDev = new double[numPoints]; // calculated Idev desmeared
    SASinstrumentType instrument = (SASinstrumentType) entry.getSASinstrument();
    SASdetectorType detector = (SASdetectorType) instrument.getSASdetector().get( ←
        detectorIndex);
    slit_length = detector.getSlitLength().getValue();
}

/**
 * @param (String) pkg Java package containing XML Schema bound to Java data structures
 * @param (String) xmlFile XML file to be opened
 * @return (Object) root object of Java data structure from XML file
 * @throws JAXBException
 */
@SuppressWarnings("unchecked")
private static Object loadXML(String pkg, String xmlFile) throws JAXBException {
    // use the $(pkg) schema that is bound to a Java structure
    JAXBContext jc = JAXBContext.newInstance(pkg);
    Unmarshaller unmarshaller = jc.createUnmarshaller();
    // open the XML file into a Java data structure
    Object obj = (Object) ((JAXBElement<Object>) unmarshaller
        .unmarshal(new File(xmlFile))).getValue();
    return obj;
}

/**
 * @param dt (DesmearingType) Desmearing properties
 * @param sasRoot (SASrootType) SAS data from XML file
 */
public void inputReporter()
{
    System.out.println("dataFile:\t" + dt.getDataFile().trim());
    System.out.printf("dataset selected:\t/SASroot/SASentry[%d]/SASdata[%d]\n",
        dt.getEntryIndex(), dt.getDataIndex());
    System.out.printf("detector selected:\t/SASroot/SASentry[%d]/SASinstrument/SASdetector ←
        [%d]\n",
        dt.getEntryIndex(), dt.getDataIndex(), dt.getDetectorIndex());
    System.out.println("extrapolation_form:\t" + dt.getExtrapolationForm().trim());
    System.out.println("x_start_extrapolation_evaluation:\t" + dt. ←
        getXStartExtrapolationEvaluation().getValue());
    System.out.println("x_start_extrapolation_evaluation unit:\t" + dt. ←
        getXStartExtrapolationEvaluation().getUnit());
    System.out.println("iterations:\t" + dt.getIterations());
    System.out.println("iterative_weight_method:\t" + dt.getIterativeWeightMethod().trim()) ←
        ;

    System.out.println("#-----");

    int numEntries = sasRoot.getSASentry().size();
    System.out.println("SASentry elements: " + numEntries);
    for( int i = 0; i < numEntries; i++ ) {
        System.out.println("SASentry");
        SASentryType entry = sasRoot.getSASentry().get(i);
        System.out.printf("Title:\t%s\n", entry.getTitle());
        List<SASentryType.Run> runs = entry.getRun();
        System.out.printf("#Runs:\t%d\n", runs.size());
        for( int j = 0; j < runs.size(); j++ ) {
            Run run = (Run) runs.get(j);
            System.out.printf("Run@name:\t%s\n", run.getName());
        }
    }
}
```

```
        System.out.printf("Run:\t%s\n", run.getValue());
    }
    List<SASdataType> datasets = entry.getSASdata();
    System.out.printf("#SASdata:\t%d\n", datasets.size());
    for( int j = 0; j < datasets.size(); j++ ) {
        SASdataType sdt = (SASdataType) datasets.get(j);
        System.out.printf("SASdata@name:\t%s\n", sdt.getName());
        System.out.printf("#points:\t%d\n", sdt.getIdata().size());
    }
    List<SASdetectorType> detectors = entry.getSASinstrument().getSASdetector();
    System.out.printf("#SASdetector:\t%d\n", detectors.size());
    for( int j = 0; j < detectors.size(); j++ ) {
        SASdetectorType det = (SASdetectorType) detectors.get(j);
        System.out.printf("SASdata@name:\t%s\n", det.getName());
        try {
            System.out.printf("SDD:\t%g\t(%)s\n", det.getSDD()
                .getValue(), det.getSDD().getUnit());
        } catch (Exception e) {
            System.out.println("SDD:\tundefined");
        }
        try {
            System.out.printf("slit_length:\t%g\t(%)s\n", det
                .getSlitLength().getValue(), det.getSlitLength()
                .getUnit());
        } catch (Exception e) {
            System.out.println("slit_length:\tundefined");
        }
    }
    System.out.println();
}

/**
 * @param args
 */
public static void main(String[] args) {
    // load test desmearing properties and data into memory
    GetSASdata sas = new GetSASdata("test.xml");
    sas.inputReporter();
    System.out.println("the end.");
}

/**
 * @return the sasRoot
 */
public SASrootType getSasRoot() {
    return sasRoot;
}

/**
 * @param sasRoot the sasRoot to set
 */
public void setSasRoot(SASrootType sasRoot) {
    GetDesmearingInfo.sasRoot = sasRoot;
}

/**
 * @return the qsas
 */
public double[] getQsas() {
    return Qsas;
}
```

```
}

/**
 * @param qsas the qsas to set
 */
public void setQsas(double[] qsas) {
    Qsas = qsas;
}

/**
 * @return the isas
 */
public double[] getIsas() {
    return Isas;
}

/**
 * @param isas the isas to set
 */
public void setIsas(double[] isas) {
    Isas = isas;
}

/**
 * @return the idev
 */
public double[] getIddev() {
    return Idev;
}

/**
 * @param idev the idev to set
 */
public void setIddev(double[] idev) {
    Idev = idev;
}

/**
 * @return the ismr
 */
public double[] getIsmr() {
    return Ismr;
}

/**
 * @param ismr the ismr to set
 */
public void setIsmr(double[] ismr) {
    Ismr = ismr;
}

/**
 * @return the idsm
 */
public double[] getIdsm() {
    return Idsm;
}

/**
 * @param idsm the idsm to set
 */
public void setIdsm(double[] idsm) {
```

```
    Idsm = idsm;
}

/**
 * @return the idsmDev
 */
public double[] getIdsmDev() {
    return IdsmDev;
}

/**
 * @param idsmDev the idsmDev to set
 */
public void setIdsmDev(double[] idsmDev) {
    IdsmDev = idsmDev;
}

/**
 * @return
 */
public double getSlitLength() {
    return slit_length;
}
}
```

### 3.3.3.2 test.xml

Ok, better to use SVN/TRAC for these files. This example will be improved but it proves the point.

```
<?xml version="1.0"?>
<?xml-stylesheet type="text/xsl" href="example.xsl" ?>
<SASroot version="1.0"
  xmlns="cansas1d/1.0"
  xmlns:xsi="http://www.w3.org/2001/XMLSchema-instance"
  xsi:schemaLocation="cansas1d/1.0 http://svn.smallangles.net/svn/canSAS/ldwg/trunk/ ↵
    cansas1d.xsd"
>
<SASentry>
  <Title>standard test dataset for Lake desmearing routine</Title>
  <Run>Run</Run>
  <SASdata name="slit-smeared">
    <Idata><Q unit="1/A">0.000371484</Q><I unit="1/cm">211554</I><Idev unit="1/cm ↵
      ">1874.86</Idev></Idata>
    <Idata><Q unit="1/A">0.000386255</Q><I unit="1/cm">201603</I><Idev unit="1/cm ↵
      ">1721.35</Idev></Idata>
    <Idata><Q unit="1/A">0.000392446</Q><I unit="1/cm">193423</I><Idev unit="1/cm ↵
      ">4250.66</Idev></Idata>
    <Idata><Q unit="1/A">0.000400937</Q><I unit="1/cm">205280</I><Idev unit="1/cm ↵
      ">1563.25</Idev></Idata>
    <Idata><Q unit="1/A">0.000415708</Q><I unit="1/cm">198569</I><Idev unit="1/cm ↵
      ">1446.58</Idev></Idata>
    <Idata><Q unit="1/A">0.000430391</Q><I unit="1/cm">198201</I><Idev unit="1/cm ↵
      ">1334.48</Idev></Idata>
    <Idata><Q unit="1/A">0.000445162</Q><I unit="1/cm">191430</I><Idev unit="1/cm ↵
      ">624.224</Idev></Idata>
    <Idata><Q unit="1/A">0.000451353</Q><I unit="1/cm">188171</I><Idev unit="1/cm ↵
      ">605.955</Idev></Idata>
    <Idata><Q unit="1/A">0.000459932</Q><I unit="1/cm">192450</I><Idev unit="1/cm ↵
      ">592.662</Idev></Idata>
```

```
<Idata><Q unit="1/A">0.000474703</Q><I unit="1/cm">186589</I><Idev unit="1/cm ↔  
">566.562</Idev></Idata>  
<Idata><Q unit="1/A">0.000489386</Q><I unit="1/cm">184247</I><Idev unit="1/cm ↔  
">541.442</Idev></Idata>  
<Idata><Q unit="1/A">0.000504157</Q><I unit="1/cm">179316</I><Idev unit="1/cm ↔  
">519.719</Idev></Idata>  
<Idata><Q unit="1/A">0.000510348</Q><I unit="1/cm">172441</I><Idev unit="1/cm ↔  
">505.368</Idev></Idata>  
<Idata><Q unit="1/A">0.000518928</Q><I unit="1/cm">175473</I><Idev unit="1/cm ↔  
">497.727</Idev></Idata>  
<Idata><Q unit="1/A">0.000533699</Q><I unit="1/cm">171012</I><Idev unit="1/cm ↔  
">479.382</Idev></Idata>  
<Idata><Q unit="1/A">0.000548381</Q><I unit="1/cm">167081</I><Idev unit="1/cm ↔  
">461.221</Idev></Idata>  
<Idata><Q unit="1/A">0.000563063</Q><I unit="1/cm">162303</I><Idev unit="1/cm ↔  
">446.693</Idev></Idata>  
<Idata><Q unit="1/A">0.000569255</Q><I unit="1/cm">158623</I><Idev unit="1/cm ↔  
">439.165</Idev></Idata>  
<Idata><Q unit="1/A">0.000577834</Q><I unit="1/cm">160015</I><Idev unit="1/cm ↔  
">434.267</Idev></Idata>  
<Idata><Q unit="1/A">0.000592605</Q><I unit="1/cm">155494</I><Idev unit="1/cm ↔  
">421.757</Idev></Idata>  
<Idata><Q unit="1/A">0.000607376</Q><I unit="1/cm">151073</I><Idev unit="1/cm ↔  
">407.866</Idev></Idata>  
<Idata><Q unit="1/A">0.000622059</Q><I unit="1/cm">146555</I><Idev unit="1/cm ↔  
">396.055</Idev></Idata>  
<Idata><Q unit="1/A">0.00062825</Q><I unit="1/cm">143885</I><Idev unit="1/cm ↔  
">390.636</Idev></Idata>  
<Idata><Q unit="1/A">0.00063683</Q><I unit="1/cm">143034</I><Idev unit="1/cm ↔  
">384.251</Idev></Idata>  
<Idata><Q unit="1/A">0.0006516</Q><I unit="1/cm">139041</I><Idev unit="1/cm ↔  
">373.826</Idev></Idata>  
<Idata><Q unit="1/A">0.000666371</Q><I unit="1/cm">136947</I><Idev unit="1/cm ↔  
">365.092</Idev></Idata>  
<Idata><Q unit="1/A">0.000681054</Q><I unit="1/cm">134324</I><Idev unit="1/cm ↔  
">357.809</Idev></Idata>  
<Idata><Q unit="1/A">0.000687245</Q><I unit="1/cm">131392</I><Idev unit="1/cm ↔  
">352.914</Idev></Idata>  
<Idata><Q unit="1/A">0.000695825</Q><I unit="1/cm">131867</I><Idev unit="1/cm ↔  
">350.159</Idev></Idata>  
<Idata><Q unit="1/A">0.000710507</Q><I unit="1/cm">128250</I><Idev unit="1/cm ↔  
">342.395</Idev></Idata>  
<Idata><Q unit="1/A">0.000725278</Q><I unit="1/cm">125404</I><Idev unit="1/cm ↔  
">334.673</Idev></Idata>  
<Idata><Q unit="1/A">0.000739961</Q><I unit="1/cm">122355</I><Idev unit="1/cm ↔  
">327.753</Idev></Idata>  
<Idata><Q unit="1/A">0.000746152</Q><I unit="1/cm">119544</I><Idev unit="1/cm ↔  
">323.132</Idev></Idata>  
<Idata><Q unit="1/A">0.000754732</Q><I unit="1/cm">118748</I><Idev unit="1/cm ↔  
">319.076</Idev></Idata>  
<Idata><Q unit="1/A">0.000769502</Q><I unit="1/cm">115545</I><Idev unit="1/cm ↔  
">311.966</Idev></Idata>  
<Idata><Q unit="1/A">0.000784273</Q><I unit="1/cm">113124</I><Idev unit="1/cm ↔  
">305.553</Idev></Idata>  
<Idata><Q unit="1/A">0.000798956</Q><I unit="1/cm">110665</I><Idev unit="1/cm ↔  
">299.932</Idev></Idata>  
<Idata><Q unit="1/A">0.000805147</Q><I unit="1/cm">109629</I><Idev unit="1/cm ↔  
">298.178</Idev></Idata>  
<Idata><Q unit="1/A">0.000813727</Q><I unit="1/cm">108497</I><Idev unit="1/cm ↔  
">294.437</Idev></Idata>  
<Idata><Q unit="1/A">0.000828498</Q><I unit="1/cm">106067</I><Idev unit="1/cm ↔  
">289.296</Idev></Idata>
```



```
<Idata><Q unit="1/A">0.000843269</Q><I unit="1/cm">103730</I><Idev unit="1/cm ↔  
">283.712</Idev></Idata>  
<Idata><Q unit="1/A">0.000857951</Q><I unit="1/cm">101055</I><Idev unit="1/cm ↔  
">278.13</Idev></Idata>  
<Idata><Q unit="1/A">0.000864142</Q><I unit="1/cm">100263</I><Idev unit="1/cm ↔  
">276.837</Idev></Idata>  
<Idata><Q unit="1/A">0.000872633</Q><I unit="1/cm">98975.8</I><Idev unit="1/cm ↔  
">272.922</Idev></Idata>  
<Idata><Q unit="1/A">0.000887404</Q><I unit="1/cm">96617</I><Idev unit="1/cm ↔  
">267.848</Idev></Idata>  
<Idata><Q unit="1/A">0.000902175</Q><I unit="1/cm">94721.9</I><Idev unit="1/cm ↔  
">263.148</Idev></Idata>  
<Idata><Q unit="1/A">0.000916946</Q><I unit="1/cm">92784.7</I><Idev unit="1/cm ↔  
">258.983</Idev></Idata>  
<Idata><Q unit="1/A">0.000923138</Q><I unit="1/cm">92237.6</I><Idev unit="1/cm ↔  
">258.352</Idev></Idata>  
<Idata><Q unit="1/A">0.000931629</Q><I unit="1/cm">91165.3</I><Idev unit="1/cm ↔  
">255.102</Idev></Idata>  
<Idata><Q unit="1/A">0.0009464</Q><I unit="1/cm">89261.8</I><Idev unit="1/cm ↔  
">251.208</Idev></Idata>  
<Idata><Q unit="1/A">0.00096117</Q><I unit="1/cm">87432.7</I><Idev unit="1/cm ↔  
">247.126</Idev></Idata>  
<Idata><Q unit="1/A">0.000975941</Q><I unit="1/cm">85430.8</I><Idev unit="1/cm ↔  
">243.068</Idev></Idata>  
<Idata><Q unit="1/A">0.000982133</Q><I unit="1/cm">84676.8</I><Idev unit="1/cm ↔  
">242.072</Idev></Idata>  
<Idata><Q unit="1/A">0.000990624</Q><I unit="1/cm">83582.8</I><Idev unit="1/cm ↔  
">238.943</Idev></Idata>  
<Idata><Q unit="1/A">0.00100539</Q><I unit="1/cm">81637.9</I><Idev unit="1/cm ↔  
">234.966</Idev></Idata>  
<Idata><Q unit="1/A">0.00102008</Q><I unit="1/cm">80007.9</I><Idev unit="1/cm ↔  
">231.251</Idev></Idata>  
<Idata><Q unit="1/A">0.00103485</Q><I unit="1/cm">78480.9</I><Idev unit="1/cm ↔  
">228.051</Idev></Idata>  
<Idata><Q unit="1/A">0.00104104</Q><I unit="1/cm">78020.3</I><Idev unit="1/cm ↔  
">227.629</Idev></Idata>  
<Idata><Q unit="1/A">0.00104953</Q><I unit="1/cm">76980.6</I><Idev unit="1/cm ↔  
">225.019</Idev></Idata>  
<Idata><Q unit="1/A">0.0010643</Q><I unit="1/cm">75445.5</I><Idev unit="1/cm ↔  
">221.885</Idev></Idata>  
<Idata><Q unit="1/A">0.00107907</Q><I unit="1/cm">73855.7</I><Idev unit="1/cm ↔  
">218.547</Idev></Idata>  
<Idata><Q unit="1/A">0.00109384</Q><I unit="1/cm">72137.1</I><Idev unit="1/cm ↔  
">215.031</Idev></Idata>  
<Idata><Q unit="1/A">0.00110003</Q><I unit="1/cm">71587.7</I><Idev unit="1/cm ↔  
">214.305</Idev></Idata>  
<Idata><Q unit="1/A">0.00110853</Q><I unit="1/cm">70608.1</I><Idev unit="1/cm ↔  
">211.782</Idev></Idata>  
<Idata><Q unit="1/A">0.0011233</Q><I unit="1/cm">69061.3</I><Idev unit="1/cm ↔  
">208.548</Idev></Idata>  
<Idata><Q unit="1/A">0.00113807</Q><I unit="1/cm">67611.7</I><Idev unit="1/cm ↔  
">205.543</Idev></Idata>  
<Idata><Q unit="1/A">0.00115284</Q><I unit="1/cm">66343.4</I><Idev unit="1/cm ↔  
">202.81</Idev></Idata>  
<Idata><Q unit="1/A">0.00115903</Q><I unit="1/cm">65983.2</I><Idev unit="1/cm ↔  
">202.392</Idev></Idata>  
<Idata><Q unit="1/A">0.00116752</Q><I unit="1/cm">65188.2</I><Idev unit="1/cm ↔  
">200.268</Idev></Idata>  
<Idata><Q unit="1/A">0.0011822</Q><I unit="1/cm">64072.9</I><Idev unit="1/cm ↔  
">197.834</Idev></Idata>  
<Idata><Q unit="1/A">0.00119697</Q><I unit="1/cm">62740.3</I><Idev unit="1/cm ↔  
">195.054</Idev></Idata>
```

```
<Idata><Q unit="1/A">0.00121175</Q><I unit="1/cm">61380.9</I><Idev unit="1/cm ↔  
">192.227</Idev></Idata>  
<Idata><Q unit="1/A">0.00121794</Q><I unit="1/cm">60799.2</I><Idev unit="1/cm ↔  
">191.38</Idev></Idata>  
<Idata><Q unit="1/A">0.00122652</Q><I unit="1/cm">60016.7</I><Idev unit="1/cm ↔  
">189.35</Idev></Idata>  
<Idata><Q unit="1/A">0.0012412</Q><I unit="1/cm">58722.4</I><Idev unit="1/cm ↔  
">186.613</Idev></Idata>  
<Idata><Q unit="1/A">0.00125597</Q><I unit="1/cm">57573.3</I><Idev unit="1/cm ↔  
">184.131</Idev></Idata>  
<Idata><Q unit="1/A">0.00127693</Q><I unit="1/cm">56176.4</I><Idev unit="1/cm ↔  
">181.407</Idev></Idata>  
<Idata><Q unit="1/A">0.00133584</Q><I unit="1/cm">52010.5</I><Idev unit="1/cm ↔  
">172.352</Idev></Idata>  
<Idata><Q unit="1/A">0.00139483</Q><I unit="1/cm">48213.1</I><Idev unit="1/cm ↔  
">163.997</Idev></Idata>  
<Idata><Q unit="1/A">0.00145383</Q><I unit="1/cm">44843.2</I><Idev unit="1/cm ↔  
">156.44</Idev></Idata>  
<Idata><Q unit="1/A">0.00151274</Q><I unit="1/cm">41576.4</I><Idev unit="1/cm ↔  
">149.07</Idev></Idata>  
<Idata><Q unit="1/A">0.00157173</Q><I unit="1/cm">38658.5</I><Idev unit="1/cm ↔  
">142.343</Idev></Idata>  
<Idata><Q unit="1/A">0.00160092</Q><I unit="1/cm">37457.7</I><Idev unit="1/cm ↔  
">139.705</Idev></Idata>  
<Idata><Q unit="1/A">0.00163073</Q><I unit="1/cm">36024.4</I><Idev unit="1/cm ↔  
">136.191</Idev></Idata>  
<Idata><Q unit="1/A">0.00168963</Q><I unit="1/cm">33573.8</I><Idev unit="1/cm ↔  
">130.321</Idev></Idata>  
<Idata><Q unit="1/A">0.00174863</Q><I unit="1/cm">31319.9</I><Idev unit="1/cm ↔  
">124.828</Idev></Idata>  
<Idata><Q unit="1/A">0.00180753</Q><I unit="1/cm">29339.2</I><Idev unit="1/cm ↔  
">119.87</Idev></Idata>  
<Idata><Q unit="1/A">0.00183672</Q><I unit="1/cm">28494.9</I><Idev unit="1/cm ↔  
">117.853</Idev></Idata>  
<Idata><Q unit="1/A">0.00186653</Q><I unit="1/cm">27449.7</I><Idev unit="1/cm ↔  
">115.05</Idev></Idata>  
<Idata><Q unit="1/A">0.00192553</Q><I unit="1/cm">25747.7</I><Idev unit="1/cm ↔  
">110.669</Idev></Idata>  
<Idata><Q unit="1/A">0.00198443</Q><I unit="1/cm">24113.1</I><Idev unit="1/cm ↔  
">106.346</Idev></Idata>  
<Idata><Q unit="1/A">0.00204352</Q><I unit="1/cm">22592.2</I><Idev unit="1/cm ↔  
">102.201</Idev></Idata>  
<Idata><Q unit="1/A">0.0020727</Q><I unit="1/cm">22004.4</I><Idev unit="1/cm ↔  
">100.667</Idev></Idata>  
<Idata><Q unit="1/A">0.00210242</Q><I unit="1/cm">21173.5</I><Idev unit="1/cm ↔  
">98.3129</Idev></Idata>  
<Idata><Q unit="1/A">0.00216142</Q><I unit="1/cm">19890</I><Idev unit="1/cm ↔  
">90.6563</Idev></Idata>  
<Idata><Q unit="1/A">0.00222041</Q><I unit="1/cm">18698.5</I><Idev unit="1/cm ↔  
">85.5865</Idev></Idata>  
<Idata><Q unit="1/A">0.00227932</Q><I unit="1/cm">17602</I><Idev unit="1/cm ↔  
">82.5591</Idev></Idata>  
<Idata><Q unit="1/A">0.00230851</Q><I unit="1/cm">17224</I><Idev unit="1/cm ↔  
">83.2993</Idev></Idata>  
<Idata><Q unit="1/A">0.00233831</Q><I unit="1/cm">16587.4</I><Idev unit="1/cm ↔  
">83.7386</Idev></Idata>  
<Idata><Q unit="1/A">0.00239722</Q><I unit="1/cm">15654</I><Idev unit="1/cm ↔  
">77.0515</Idev></Idata>  
<Idata><Q unit="1/A">0.00245622</Q><I unit="1/cm">14772.5</I><Idev unit="1/cm ↔  
">74.475</Idev></Idata>  
<Idata><Q unit="1/A">0.00251521</Q><I unit="1/cm">13950.9</I><Idev unit="1/cm ↔  
">72.0371</Idev></Idata>
```

```
<Idata><Q unit="1/A">0.0025444</Q><I unit="1/cm">13671.6</I><Idev unit="1/cm ↔  
">71.2492</Idev></Idata>  
<Idata><Q unit="1/A">0.00257412</Q><I unit="1/cm">13179.3</I><Idev unit="1/cm ↔  
">69.6967</Idev></Idata>  
<Idata><Q unit="1/A">0.00263311</Q><I unit="1/cm">12441.4</I><Idev unit="1/cm ↔  
">67.4163</Idev></Idata>  
<Idata><Q unit="1/A">0.00269211</Q><I unit="1/cm">11793.6</I><Idev unit="1/cm ↔  
">65.3736</Idev></Idata>  
<Idata><Q unit="1/A">0.0027802</Q><I unit="1/cm">11030.1</I><Idev unit="1/cm ↔  
">62.915</Idev></Idata>  
<Idata><Q unit="1/A">0.0030161</Q><I unit="1/cm">8974.34</I><Idev unit="1/cm ↔  
">55.8783</Idev></Idata>  
<Idata><Q unit="1/A">0.00325199</Q><I unit="1/cm">7394.04</I><Idev unit="1/cm ↔  
">50.0284</Idev></Idata>  
<Idata><Q unit="1/A">0.00348788</Q><I unit="1/cm">6118.83</I><Idev unit="1/cm ↔  
">44.9183</Idev></Idata>  
<Idata><Q unit="1/A">0.00372368</Q><I unit="1/cm">5113.03</I><Idev unit="1/cm ↔  
">40.5794</Idev></Idata>  
<Idata><Q unit="1/A">0.00395958</Q><I unit="1/cm">4318.86</I><Idev unit="1/cm ↔  
">7.43972</Idev></Idata>  
<Idata><Q unit="1/A">0.00419538</Q><I unit="1/cm">3658.3</I><Idev unit="1/cm ↔  
">6.64065</Idev></Idata>  
<Idata><Q unit="1/A">0.00443127</Q><I unit="1/cm">3139.32</I><Idev unit="1/cm ↔  
">5.99763</Idev></Idata>  
<Idata><Q unit="1/A">0.00466716</Q><I unit="1/cm">2713.56</I><Idev unit="1/cm ↔  
">5.45451</Idev></Idata>  
<Idata><Q unit="1/A">0.00490297</Q><I unit="1/cm">2365</I><Idev unit="1/cm">5.00216</ ↔  
Idev></Idata>  
<Idata><Q unit="1/A">0.00513886</Q><I unit="1/cm">2070.02</I><Idev unit="1/cm ↔  
">4.60551</Idev></Idata>  
<Idata><Q unit="1/A">0.00537466</Q><I unit="1/cm">1830.58</I><Idev unit="1/cm ↔  
">4.28091</Idev></Idata>  
<Idata><Q unit="1/A">0.00561065</Q><I unit="1/cm">1625.29</I><Idev unit="1/cm ↔  
">3.99215</Idev></Idata>  
<Idata><Q unit="1/A">0.00584645</Q><I unit="1/cm">1453.64</I><Idev unit="1/cm ↔  
">3.74405</Idev></Idata>  
<Idata><Q unit="1/A">0.00608234</Q><I unit="1/cm">1305.37</I><Idev unit="1/cm ↔  
">3.52884</Idev></Idata>  
<Idata><Q unit="1/A">0.00631815</Q><I unit="1/cm">1174.14</I><Idev unit="1/cm ↔  
">3.32995</Idev></Idata>  
<Idata><Q unit="1/A">0.00655404</Q><I unit="1/cm">1064.44</I><Idev unit="1/cm ↔  
">3.16037</Idev></Idata>  
<Idata><Q unit="1/A">0.00678993</Q><I unit="1/cm">967.427</I><Idev unit="1/cm ↔  
">3.01166</Idev></Idata>  
<Idata><Q unit="1/A">0.00702582</Q><I unit="1/cm">878.938</I><Idev unit="1/cm ↔  
">2.86918</Idev></Idata>  
<Idata><Q unit="1/A">0.00726163</Q><I unit="1/cm">803.771</I><Idev unit="1/cm ↔  
">2.74895</Idev></Idata>  
<Idata><Q unit="1/A">0.00749743</Q><I unit="1/cm">734.672</I><Idev unit="1/cm ↔  
">2.63429</Idev></Idata>  
<Idata><Q unit="1/A">0.00773332</Q><I unit="1/cm">677.548</I><Idev unit="1/cm ↔  
">2.53748</Idev></Idata>  
<Idata><Q unit="1/A">0.00796921</Q><I unit="1/cm">626.597</I><Idev unit="1/cm ↔  
">2.44863</Idev></Idata>  
<Idata><Q unit="1/A">0.00820511</Q><I unit="1/cm">579.635</I><Idev unit="1/cm ↔  
">2.36697</Idev></Idata>  
<Idata><Q unit="1/A">0.00844091</Q><I unit="1/cm">536.91</I><Idev unit="1/cm ↔  
">2.29218</Idev></Idata>  
<Idata><Q unit="1/A">0.0086768</Q><I unit="1/cm">500.271</I><Idev unit="1/cm ↔  
">2.22717</Idev></Idata>  
<Idata><Q unit="1/A">0.00891261</Q><I unit="1/cm">464.949</I><Idev unit="1/cm ↔  
">2.1601</Idev></Idata>
```

```
<Idata><Q unit="1/A">0.00914859</Q><I unit="1/cm">430.058</I><Idev unit="1/cm ↔  
">2.09694</Idev></Idata>  
<Idata><Q unit="1/A">0.00938439</Q><I unit="1/cm">399.444</I><Idev unit="1/cm ↔  
">2.03912</Idev></Idata>  
<Idata><Q unit="1/A">0.00962028</Q><I unit="1/cm">375.619</I><Idev unit="1/cm ↔  
">1.9894</Idev></Idata>  
<Idata><Q unit="1/A">0.00985609</Q><I unit="1/cm">348.697</I><Idev unit="1/cm ↔  
">1.93717</Idev></Idata>  
<Idata><Q unit="1/A">0.010092</Q><I unit="1/cm">328.261</I><Idev unit="1/cm">1.8957</ ↔  
Idev></Idata>  
<Idata><Q unit="1/A">0.0103279</Q><I unit="1/cm">309.772</I><Idev unit="1/cm ↔  
">1.85859</Idev></Idata>  
<Idata><Q unit="1/A">0.0104458</Q><I unit="1/cm">302.426</I><Idev unit="1/cm ↔  
">1.84512</Idev></Idata>  
<Idata><Q unit="1/A">0.0105638</Q><I unit="1/cm">292.196</I><Idev unit="1/cm ↔  
">1.82164</Idev></Idata>  
<Idata><Q unit="1/A">0.0107996</Q><I unit="1/cm">272.191</I><Idev unit="1/cm ↔  
">1.78264</Idev></Idata>  
<Idata><Q unit="1/A">0.0110354</Q><I unit="1/cm">261.342</I><Idev unit="1/cm ↔  
">1.75923</Idev></Idata>  
<Idata><Q unit="1/A">0.0112713</Q><I unit="1/cm">245.999</I><Idev unit="1/cm ↔  
">1.72526</Idev></Idata>  
<Idata><Q unit="1/A">0.0113893</Q><I unit="1/cm">241.462</I><Idev unit="1/cm ↔  
">1.7191</Idev></Idata>  
<Idata><Q unit="1/A">0.0115072</Q><I unit="1/cm">233.87</I><Idev unit="1/cm ↔  
">1.70363</Idev></Idata>  
<Idata><Q unit="1/A">0.011743</Q><I unit="1/cm">222.267</I><Idev unit="1/cm ↔  
">1.67677</Idev></Idata>  
<Idata><Q unit="1/A">0.0119789</Q><I unit="1/cm">213.558</I><Idev unit="1/cm ↔  
">1.65675</Idev></Idata>  
<Idata><Q unit="1/A">0.0122147</Q><I unit="1/cm">202.04</I><Idev unit="1/cm ↔  
">1.63174</Idev></Idata>  
<Idata><Q unit="1/A">0.0123326</Q><I unit="1/cm">196.924</I><Idev unit="1/cm ↔  
">1.62249</Idev></Idata>  
<Idata><Q unit="1/A">0.0124505</Q><I unit="1/cm">192.761</I><Idev unit="1/cm ↔  
">1.6123</Idev></Idata>  
<Idata><Q unit="1/A">0.0126865</Q><I unit="1/cm">185.892</I><Idev unit="1/cm ↔  
">1.59555</Idev></Idata>  
<Idata><Q unit="1/A">0.0129223</Q><I unit="1/cm">176.483</I><Idev unit="1/cm ↔  
">1.57692</Idev></Idata>  
<Idata><Q unit="1/A">0.0131582</Q><I unit="1/cm">171.938</I><Idev unit="1/cm ↔  
">1.56441</Idev></Idata>  
<Idata><Q unit="1/A">0.0132761</Q><I unit="1/cm">167.199</I><Idev unit="1/cm ↔  
">1.55535</Idev></Idata>  
<Idata><Q unit="1/A">0.0142196</Q><I unit="1/cm">144.551</I><Idev unit="1/cm ↔  
">1.50508</Idev></Idata>  
<Idata><Q unit="1/A">0.0151631</Q><I unit="1/cm">126.664</I><Idev unit="1/cm ↔  
">1.46166</Idev></Idata>  
<Idata><Q unit="1/A">0.0161065</Q><I unit="1/cm">112.514</I><Idev unit="1/cm ↔  
">1.42482</Idev></Idata>  
<Idata><Q unit="1/A">0.0170499</Q><I unit="1/cm">98.3946</I><Idev unit="1/cm ↔  
">1.39316</Idev></Idata>  
<Idata><Q unit="1/A">0.0179933</Q><I unit="1/cm">90.2142</I><Idev unit="1/cm ↔  
">1.37318</Idev></Idata>  
<Idata><Q unit="1/A">0.0189368</Q><I unit="1/cm">82.8805</I><Idev unit="1/cm ↔  
">1.35248</Idev></Idata>  
<Idata><Q unit="1/A">0.0198803</Q><I unit="1/cm">75.2953</I><Idev unit="1/cm ↔  
">1.33271</Idev></Idata>  
<Idata><Q unit="1/A">0.0208238</Q><I unit="1/cm">71.261</I><Idev unit="1/cm">1.3229</ ↔  
Idev></Idata>  
<Idata><Q unit="1/A">0.0217672</Q><I unit="1/cm">64.006</I><Idev unit="1/cm ↔  
">1.30757</Idev></Idata>
```

```
<Idata><Q unit="1/A">0.0227106</Q><I unit="1/cm">61.7542</I><Idev unit="1/cm ↔  
">1.30381</Idev></Idata>  
<Idata><Q unit="1/A">0.023654</Q><I unit="1/cm">61.7168</I><Idev unit="1/cm ↔  
">1.30093</Idev></Idata>  
<Idata><Q unit="1/A">0.0245975</Q><I unit="1/cm">57.8197</I><Idev unit="1/cm ↔  
">1.28888</Idev></Idata>  
<Idata><Q unit="1/A">0.025541</Q><I unit="1/cm">54.3294</I><Idev unit="1/cm ↔  
">1.27828</Idev></Idata>  
<Idata><Q unit="1/A">0.0264845</Q><I unit="1/cm">53.7715</I><Idev unit="1/cm ↔  
">1.27858</Idev></Idata>  
<Idata><Q unit="1/A">0.0274279</Q><I unit="1/cm">51.3464</I><Idev unit="1/cm ↔  
">1.27158</Idev></Idata>  
<Idata><Q unit="1/A">0.0283713</Q><I unit="1/cm">50.7223</I><Idev unit="1/cm ↔  
">1.27033</Idev></Idata>  
<Idata><Q unit="1/A">0.0293147</Q><I unit="1/cm">48.6453</I><Idev unit="1/cm ↔  
">1.27083</Idev></Idata>  
<Idata><Q unit="1/A">0.0302582</Q><I unit="1/cm">46.8375</I><Idev unit="1/cm ↔  
">1.26587</Idev></Idata>  
<Idata><Q unit="1/A">0.0312017</Q><I unit="1/cm">47.594</I><Idev unit="1/cm ↔  
">1.26671</Idev></Idata>  
<Idata><Q unit="1/A">0.0321451</Q><I unit="1/cm">44.9242</I><Idev unit="1/cm ↔  
">1.25974</Idev></Idata>  
<Idata><Q unit="1/A">0.0330885</Q><I unit="1/cm">42.9397</I><Idev unit="1/cm ↔  
">1.25796</Idev></Idata>  
<Idata><Q unit="1/A">0.034032</Q><I unit="1/cm">44.3886</I><Idev unit="1/cm ↔  
">1.25858</Idev></Idata>  
<Idata><Q unit="1/A">0.0349755</Q><I unit="1/cm">44.6934</I><Idev unit="1/cm ↔  
">1.25971</Idev></Idata>  
<Idata><Q unit="1/A">0.0359189</Q><I unit="1/cm">44.6929</I><Idev unit="1/cm ↔  
">1.26103</Idev></Idata>  
<Idata><Q unit="1/A">0.0368623</Q><I unit="1/cm">43.0895</I><Idev unit="1/cm ↔  
">1.25534</Idev></Idata>  
<Idata><Q unit="1/A">0.0378057</Q><I unit="1/cm">43.2662</I><Idev unit="1/cm ↔  
">1.25507</Idev></Idata>  
<Idata><Q unit="1/A">0.0387492</Q><I unit="1/cm">42.1147</I><Idev unit="1/cm ↔  
">1.25495</Idev></Idata>  
<Idata><Q unit="1/A">0.0396927</Q><I unit="1/cm">41.2501</I><Idev unit="1/cm ↔  
">1.25071</Idev></Idata>  
<Idata><Q unit="1/A">0.0406362</Q><I unit="1/cm">41.5693</I><Idev unit="1/cm ↔  
">1.25334</Idev></Idata>  
<Idata><Q unit="1/A">0.0415795</Q><I unit="1/cm">41.4826</I><Idev unit="1/cm ↔  
">1.25233</Idev></Idata>  
<Idata><Q unit="1/A">0.0425229</Q><I unit="1/cm">42.423</I><Idev unit="1/cm ↔  
">1.25764</Idev></Idata>  
<Idata><Q unit="1/A">0.0429947</Q><I unit="1/cm">40.4159</I><Idev unit="1/cm ↔  
">1.25491</Idev></Idata>  
<Idata><Q unit="1/A">0.0434664</Q><I unit="1/cm">41.3698</I><Idev unit="1/cm ↔  
">1.25613</Idev></Idata>  
<Idata><Q unit="1/A">0.0444099</Q><I unit="1/cm">39.2216</I><Idev unit="1/cm ↔  
">1.25011</Idev></Idata>  
<Idata><Q unit="1/A">0.0453533</Q><I unit="1/cm">40.7132</I><Idev unit="1/cm ↔  
">1.25594</Idev></Idata>  
<Idata><Q unit="1/A">0.0462968</Q><I unit="1/cm">40.6365</I><Idev unit="1/cm ↔  
">1.25534</Idev></Idata>  
<Idata><Q unit="1/A">0.0467685</Q><I unit="1/cm">40.1072</I><Idev unit="1/cm ↔  
">1.25637</Idev></Idata>  
<Idata><Q unit="1/A">0.0472402</Q><I unit="1/cm">39.9715</I><Idev unit="1/cm ↔  
">1.25581</Idev></Idata>  
<Idata><Q unit="1/A">0.0481837</Q><I unit="1/cm">41.5141</I><Idev unit="1/cm ↔  
">1.26048</Idev></Idata>  
<Idata><Q unit="1/A">0.0491271</Q><I unit="1/cm">39.4205</I><Idev unit="1/cm ↔  
">1.25422</Idev></Idata>
```

```
<Idata><Q unit="1/A">0.0500705</Q><I unit="1/cm">39.967</I><Idev unit="1/cm ↔  
">1.25656</Idev></Idata>  
<Idata><Q unit="1/A">0.0505422</Q><I unit="1/cm">40.7231</I><Idev unit="1/cm ↔  
">1.25899</Idev></Idata>  
<Idata><Q unit="1/A">0.051014</Q><I unit="1/cm">39.0016</I><Idev unit="1/cm ↔  
">1.25301</Idev></Idata>  
<Idata><Q unit="1/A">0.0519574</Q><I unit="1/cm">38.1899</I><Idev unit="1/cm ↔  
">1.25023</Idev></Idata>  
<Idata><Q unit="1/A">0.0529008</Q><I unit="1/cm">40.2931</I><Idev unit="1/cm ↔  
">1.2583</Idev></Idata>  
<Idata><Q unit="1/A">0.0538443</Q><I unit="1/cm">38.8024</I><Idev unit="1/cm ↔  
">1.25495</Idev></Idata>  
<Idata><Q unit="1/A">0.054316</Q><I unit="1/cm">39.2194</I><Idev unit="1/cm ↔  
">1.25647</Idev></Idata>  
<Idata><Q unit="1/A">0.0547878</Q><I unit="1/cm">37.9188</I><Idev unit="1/cm ↔  
">1.25364</Idev></Idata>  
<Idata><Q unit="1/A">0.0557312</Q><I unit="1/cm">36.8598</I><Idev unit="1/cm ↔  
">1.25236</Idev></Idata>  
<Idata><Q unit="1/A">0.0566745</Q><I unit="1/cm">38.9685</I><Idev unit="1/cm ↔  
">1.25924</Idev></Idata>  
<Idata><Q unit="1/A">0.0580898</Q><I unit="1/cm">38.4669</I><Idev unit="1/cm ↔  
">1.25824</Idev></Idata>  
<Idata><Q unit="1/A">0.0618635</Q><I unit="1/cm">39.8942</I><Idev unit="1/cm ↔  
">1.26309</Idev></Idata>  
<Idata><Q unit="1/A">0.0656372</Q><I unit="1/cm">39.609</I><Idev unit="1/cm ↔  
">1.26249</Idev></Idata>  
<Idata><Q unit="1/A">0.0694109</Q><I unit="1/cm">38.9038</I><Idev unit="1/cm ↔  
">1.2647</Idev></Idata>  
<Idata><Q unit="1/A">0.0731847</Q><I unit="1/cm">38.208</I><Idev unit="1/cm ↔  
">1.26368</Idev></Idata>  
<Idata><Q unit="1/A">0.0769583</Q><I unit="1/cm">38.1583</I><Idev unit="1/cm ↔  
">1.26555</Idev></Idata>  
<Idata><Q unit="1/A">0.0807321</Q><I unit="1/cm">39.239</I><Idev unit="1/cm ↔  
">1.26895</Idev></Idata>  
<Idata><Q unit="1/A">0.0845057</Q><I unit="1/cm">38.8689</I><Idev unit="1/cm ↔  
">1.27167</Idev></Idata>  
<Idata><Q unit="1/A">0.0882794</Q><I unit="1/cm">36.3077</I><Idev unit="1/cm ↔  
">1.26913</Idev></Idata>  
<Idata><Q unit="1/A">0.0920531</Q><I unit="1/cm">37.3417</I><Idev unit="1/cm ↔  
">1.27279</Idev></Idata>  
<Idata><Q unit="1/A">0.0958267</Q><I unit="1/cm">39.2128</I><Idev unit="1/cm ↔  
">1.27939</Idev></Idata>  
<Idata><Q unit="1/A">0.0996003</Q><I unit="1/cm">38.8772</I><Idev unit="1/cm ↔  
">1.27929</Idev></Idata>  
<Idata><Q unit="1/A">0.103374</Q><I unit="1/cm">38.2984</I><Idev unit="1/cm ↔  
">1.28285</Idev></Idata>  
<Idata><Q unit="1/A">0.107148</Q><I unit="1/cm">37.2503</I><Idev unit="1/cm ↔  
">1.28113</Idev></Idata>  
<Idata><Q unit="1/A">0.110921</Q><I unit="1/cm">38.2859</I><Idev unit="1/cm ↔  
">1.28426</Idev></Idata>  
<Idata><Q unit="1/A">0.114695</Q><I unit="1/cm">37.1071</I><Idev unit="1/cm ↔  
">1.28385</Idev></Idata>  
<Idata><Q unit="1/A">0.118468</Q><I unit="1/cm">37.2796</I><Idev unit="1/cm ↔  
">1.28837</Idev></Idata>  
<Idata><Q unit="1/A">0.122242</Q><I unit="1/cm">38.1078</I><Idev unit="1/cm ↔  
">1.29084</Idev></Idata>  
<Idata><Q unit="1/A">0.126015</Q><I unit="1/cm">37.9901</I><Idev unit="1/cm ↔  
">1.29471</Idev></Idata>  
<Idata><Q unit="1/A">0.129789</Q><I unit="1/cm">37.2332</I><Idev unit="1/cm ↔  
">1.29618</Idev></Idata>  
<Idata><Q unit="1/A">0.133562</Q><I unit="1/cm">37.525</I><Idev unit="1/cm">1.29856</ ↔  
Idev></Idata>
```

```
<Idata><Q unit="1/A">0.137336</Q><I unit="1/cm">39.7959</I><Idev unit="1/cm ↔  
">1.30577</Idev></Idata>  
<Idata><Q unit="1/A">0.141109</Q><I unit="1/cm">37.5901</I><Idev unit="1/cm ↔  
">1.30129</Idev></Idata>  
<Idata><Q unit="1/A">0.144883</Q><I unit="1/cm">37.5137</I><Idev unit="1/cm ↔  
">1.30765</Idev></Idata>  
<Idata><Q unit="1/A">0.148656</Q><I unit="1/cm">38.3692</I><Idev unit="1/cm ↔  
">1.31024</Idev></Idata>  
<Idata><Q unit="1/A">0.152429</Q><I unit="1/cm">37.9165</I><Idev unit="1/cm ↔  
">1.30974</Idev></Idata>  
<Idata><Q unit="1/A">0.156203</Q><I unit="1/cm">38.4753</I><Idev unit="1/cm ↔  
">1.31639</Idev></Idata>  
<Idata><Q unit="1/A">0.159976</Q><I unit="1/cm">38.8267</I><Idev unit="1/cm">1.3178</ ↔  
Idev></Idata>  
<Idata><Q unit="1/A">0.163749</Q><I unit="1/cm">37.9845</I><Idev unit="1/cm ↔  
">1.32281</Idev></Idata>  
<Idata><Q unit="1/A">0.167523</Q><I unit="1/cm">39.9222</I><Idev unit="1/cm ↔  
">1.32564</Idev></Idata>  
<Idata><Q unit="1/A">0.171296</Q><I unit="1/cm">41.1806</I><Idev unit="1/cm ↔  
">1.32861</Idev></Idata>  
<Idata><Q unit="1/A">0.175069</Q><I unit="1/cm">38.5425</I><Idev unit="1/cm ↔  
">1.32443</Idev></Idata>  
<Idata><Q unit="1/A">0.178843</Q><I unit="1/cm">39.2107</I><Idev unit="1/cm ↔  
">1.33117</Idev></Idata>  
<Idata><Q unit="1/A">0.182616</Q><I unit="1/cm">38.3168</I><Idev unit="1/cm ↔  
">1.33153</Idev></Idata>  
<Idata><Q unit="1/A">0.186389</Q><I unit="1/cm">40.2098</I><Idev unit="1/cm ↔  
">1.33532</Idev></Idata>  
<Idata><Q unit="1/A">0.190162</Q><I unit="1/cm">39.1407</I><Idev unit="1/cm ↔  
">1.33442</Idev></Idata>  
<Idata><Q unit="1/A">0.193935</Q><I unit="1/cm">38.3557</I><Idev unit="1/cm ↔  
">1.33652</Idev></Idata>  
<Idata><Q unit="1/A">0.197708</Q><I unit="1/cm">39.7318</I><Idev unit="1/cm ↔  
">1.34276</Idev></Idata>  
<Idata><Q unit="1/A">0.201481</Q><I unit="1/cm">36.7008</I><Idev unit="1/cm ↔  
">1.33804</Idev></Idata>  
<Idata><Q unit="1/A">0.205254</Q><I unit="1/cm">37.2223</I><Idev unit="1/cm ↔  
">1.34263</Idev></Idata>  
<Idata><Q unit="1/A">0.209027</Q><I unit="1/cm">39.6126</I><Idev unit="1/cm ↔  
">1.34766</Idev></Idata>  
<Idata><Q unit="1/A">0.2128</Q><I unit="1/cm">37.604</I><Idev unit="1/cm">1.34668</ ↔  
Idev></Idata>  
<Idata><Q unit="1/A">0.216573</Q><I unit="1/cm">39.0708</I><Idev unit="1/cm">1.3538</ ↔  
Idev></Idata>  
<Idata><Q unit="1/A">0.220346</Q><I unit="1/cm">38.2783</I><Idev unit="1/cm ↔  
">1.35074</Idev></Idata>  
<Idata><Q unit="1/A">0.224119</Q><I unit="1/cm">38.589</I><Idev unit="1/cm">1.35581</ ↔  
Idev></Idata>  
</SASdata>  
<SASsample>  
  <ID>ID</ID>  
</SASsample>  
<SASinstrument>  
  <name>calculated</name>  
  <SASsource>  
    <radiation>model</radiation>  
  </SASsource>  
  <SAScollimation/>  
  <SASdetector>  
    <name>calculated</name>  
    <slit_length unit="1/A">0.08</slit_length>  
  </SASdetector>
```

```
</SASinstrument>
<SASnote/>
</SASentry>
</SASroot>
```

## 3.4 Python binding

Specific *support* for the cansas1d/1.0 data standard in Python is being developed by *NIST/NCNR* as part of their contribution to the *DANSE* project.

Here are some terse instructions to get you started:

```
svn co http://danse.us/trac/sans/browser/trunk/DataLoader DataLoader
cd DataLoader
python setup.py install
```

The *release notes* have a list of the dependencies.

### 3.4.1 Comments

Other constructive suggestions (that predate the NIST/NCNR support) have been gathered on this page.

### 3.4.2 gnosis.xml.objectify

The GnosisUtils (<http://freshmeat.net/projects/gnosisxml/>) offer a method to read any XML file into Python data structures. This utility does not validate the XML against a specific XML Schema which can be both good (flexible, especially when XML Foreign Namespace elements are used) and not so good (XML content not guaranteed to be *by the rules*).

A quick test of this turned up an acceptable result in that it was able to read several of the canSAS test XML files, including those with foreign namespaces. And *it was very easy*. (Especially with some help from <http://www.xml.com/pub/a/2003/07/02/py-xml.html>)

Here is a quick example.

#### 3.4.2.1 installation

Here is the condensed installation (without all that output) steps. Your system may have gnosis already installed. You may also need sysAdmin privileges. You may need ...

```
cd /tmp
wget http://freshmeat.net/redirect/gnosisxml/22028/url_tgz/Gnosis_Utils-1.2.2.tar.gz
tar xzf Gnosis_Utils-1.2.2.tar.gz
cd Gnosis_Utils-1.2.2/
python setup.py install_all
```

#### 3.4.2.2 quick test in Python

Here is the Python code (without all that output) (called *python-test.py*):

```
import gnosis.xml.objectify

sasxml = gnosis.xml.objectify.XML_Objectify('bimodal-test1.xml').make_instance()
print sasxml.SASentry.Title.PCDATA
print sasxml.SASentry.Run.PCDATA
print sasxml.SASentry.SASinstrument.name.PCDATA
```



```
print sasxml.SASentry.SASdata.Idata[0].Q.unit, sasxml.SASentry.SASdata.Idata[0].I.unit
print sasxml.SASentry.SASdata.Idata[0].Q.PCDATA, sasxml.SASentry.SASdata.Idata[0].I.PCDATA ←
, sasxml.SASentry.SASdata.Idata[0].Idev.PCDATA
```

### 3.4.2.3 execution of that Python code

```
[Pete@regitte,2441,cansasldwg-regitte]$ ./python-test.py
SAS bimodal test1
1992
simulated SAS calculation
1/A 1/cm
0.0040157139 3497.473 90.72816
```

### 3.4.2.4 full session output

```
[Pete@regitte,2429,/tmp]$ cd /tmp
/tmp
[Pete@regitte,2430,/tmp]$ wget http://freshmeat.net/redir/gnosisxml/22028/url_tgz/ ←
Gnosis_Utils-1.2.2.tar.gz
--11:43:16-- http://freshmeat.net/redir/gnosisxml/22028/url_tgz/Gnosis_Utils-1.2.2.tar.gz
=> 'Gnosis_Utils-1.2.2.tar.gz'
Resolving freshmeat.net... 66.35.250.168
Connecting to freshmeat.net|66.35.250.168|:80... connected.
HTTP request sent, awaiting response... 302 Found
Location: http://www.gnosis.cx/download/Gnosis_Utils.More/Gnosis_Utils-1.2.2.tar.gz [ ←
following]
--11:43:16-- http://www.gnosis.cx/download/Gnosis_Utils.More/Gnosis_Utils-1.2.2.tar.gz
=> 'Gnosis_Utils-1.2.2.tar.gz'
Resolving www.gnosis.cx... 64.41.64.172
Connecting to www.gnosis.cx|64.41.64.172|:80... connected.
HTTP request sent, awaiting response... 200 OK
Length: 287,989 (281K) [application/x-tar]

100%[=====
287,989 --.-K/s

11:43:16 (2.47 MB/s) - 'Gnosis_Utils-1.2.2.tar.gz' saved [287989/287989]

[Pete@regitte,2431,/tmp]$ tar xzf Gnosis_Utils-1.2.2.tar.gz
[Pete@regitte,2432,/tmp]$ cd Gnosis_Utils-1.2.2/
/tmp/Gnosis_Utils-1.2.2
[Pete@regitte,2433,Gnosis_Utils-1.2.2]$ python setup.py install_all
[Pete@regitte,2434,Gnosis_Utils-1.2.2]$ cd ~/workspace/cansasldwg-regitte
[Pete@regitte,2435,cansasldwg-regitte]$ python
Python 2.5.1 (r251:54863, May 18 2007, 16:56:43)
[GCC 3.4.4 (cygming special, gdc 0.12, using dmd 0.125)] on cygwin
Type "help", "copyright", "credits" or "license" for more information.
>>> import gnosis.xml.objectify
>>> sasxml = gnosis.xml.objectify.XML_Objectify('bimodal-test1.xml').make_instance()
>>> print sasxml.SASentry.Title.PCDATA
SAS bimodal test1
>>> print sasxml.SASentry.Run.PCDATA
1992
>>> print sasxml.SASentry.SASinstrument.name.PCDATA
simulated SAS calculation
>>> print sasxml.SASentry.SASdata.Idata[0].Q.unit
1/A
>>> print sasxml.SASentry.SASdata.Idata[0].I.unit
```

```
1/cm
>>> print sasxml.SASentry.SASdata.Idata[0].Q.PCDATA, sasxml.SASentry.SASdata.Idata[0].I. ←
      PCDATA, sasxml.SASentry.SASdata.Idata[0].Idev.PCDATA
0.0040157139 3497.473 90.72816
```

#### 3.4.2.5 Conclusion: OK

This has the promise of being a useful approach to reading this format in Python. Now, how to write back out... ?

### 3.4.3 generateDS.py

*generateDS.py* (<http://www.rexx.com/~dkuhlman/>, <http://www.rexx.com/~dkuhlman/generateDS.html>) can build a binding (map the structure of the XML file directly into a Python data structure) for Python from an XML Schema. However, the cansas1d/1.0 XML schema (cansas1d.xsd) does not seem to fit the model. It seems, for now, that *generateDS-1.12a* fails on a certain *annotate* line.

```
[Pete@regitte,2402,cansas1dwg-regitte]$ python ~/generateDS-1.12a/generateDS.py -p CS1D_ -o ←
cansas1d.py -s cansas1dsubs.py cansas1d.xsd
Traceback (most recent call last):
  File "/home/Pete/generateDS-1.12a/generateDS.py", line 3997, in <module>
    main()
  File "/home/Pete/generateDS-1.12a/generateDS.py", line 3993, in main
    processIncludes, superModule=superModule)
  File "/home/Pete/generateDS-1.12a/generateDS.py", line 3909, in parseAndGenerate
    root.annotate()
AttributeError: 'NoneType' object has no attribute 'annotate'
```

#### 3.4.3.1 Conclusion: not ready yet

Either the canSAS standard (by means of the cansas1d.xsd XML Schema) is not ready or *generateDS.py* does not cover the XML Schema requirements we have at this time. Either way, this is not a viable tool to use now. (2008-05-16)

### 3.4.4 Other suggestions

- <http://www.devx.com/ibm/Article/20261>
- <http://mail.python.org/pipermail/xml-sig/2002-April/007559.html>
- <http://pywebsvcs.sourceforge.net/>

## Chapter 4

# Other matters

Various topics have been considered or presented in considering this standard. Some are described below.

### 4.1 XML Help

Listed below are various references useful in learning XML and related topics.

- **XML:** eXtensible Markup Language
    - <http://www.w3schools.com/xml/>
    - <http://www.w3.org/XML/>
    - <http://en.wikipedia.org/wiki/XML>
    - <http://www.zvon.org/xxl/XPathTutorial/General/examples.html>
  - **XSL (or XSLT):** eXtensible Stylesheet Language (Transformation)
    - <http://www.w3schools.com/xsl/>
    - <http://www.w3.org/Style/XSL/>
    - [http://en.wikipedia.org/wiki/Extensible\\_Stylesheet\\_Language](http://en.wikipedia.org/wiki/Extensible_Stylesheet_Language)
    - <http://en.wikipedia.org/wiki/XSLT>
  - **XPath:** XPath is a language for finding information in an XML document.
    - <http://www.w3schools.com/xpath/>
    - <http://www.w3.org/Style/XSL/>
    - <http://en.wikipedia.org/wiki/XPath>
  - **Schema:** An XML Schema describes the structure of an XML document.
    - <http://www.w3schools.com/schema/>
    - <http://www.w3.org/XML/Schema>
    - <http://en.wikipedia.org/wiki/XSD>
  - **XML Namespaces:** XML namespaces are used for providing uniquely named elements and attributes in an XML instance.
    - <http://www.zvon.org/xxl/NamespaceTutorial/Output>
    - [http://en.wikipedia.org/wiki/XML\\_namespaces](http://en.wikipedia.org/wiki/XML_namespaces)
    - [http://www.w3schools.com/XML/xml\\_namespaces.asp](http://www.w3schools.com/XML/xml_namespaces.asp)
-

- **XML Foreign Elements:** Inclusion of elements, at select locations, that are not defined by the cansas1d.xsd XML Schema
  - <http://books.xmlschemata.org/relaxng/relax-CHP-11-SECT-4.html>
  - <http://www.w3.org/TR/SVG/extend.html>
  - <http://www.google.com/search?q=XML+foreign+elements>